



Vogel- Johnson Agar Base w/o Tellurite (V.J. Agar)

M023

Intended Use:

Recommended for selective isolation of coagulase positive, mannitol fermenting *Staphylococcus aureus* from heavily contaminated food and clinical specimens.

Composition**

Ingredients	g / L
Tryptone	10.000
Yeast extract	5.000
Mannitol	10.000
Dipotassium hydrogen phosphate	5.000
Lithium chloride	5.000
Glycine	10.000
Phenol red	0.025
Agar	16.000
Final pH (at 25°C)	7.2±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 61.02 grams in 1000 ml purified / distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C and add 20 ml of sterile 1% Potassium Tellurite solution (FD052). Mix gently and pour into sterile Petri plates.

Principle And Interpretation

Staphylococcus aureus, a gram-positive, spherical bacterium, is a common colonizer of the human skin and mucosa. It causes skin and wound infections, urinary tract infections, pneumonia and bacteremia. It is also commonly implicated in food poisoning. It is also found as a common contaminant in pharmaceutical and cosmetics products (1). Vogel-Johnson Agar is prepared according to the formula devised by Vogel and Johnson (2) and is recommended for the microbial limit test in USP (3). Originally it was developed by Zebovitz (4), as a Tellurite Glycine Agar, a selective medium for the detection of coagulase-positive staphylococci. Vogel-Johnson modified the medium in 1960 by the addition of phenol red as a pH indicator and by increasing the quantity of mannitol (2). Selection and differentiation of coagulase-positive staphylococci on V.J. Agar is based on mannitol fermentation and tellurite reduction (5). V.J. Agar is specified in the standard methods for examination of cosmetics (1,6), pharmaceutical articles and nutritional supplements (3). In addition, the formulation complies with recommendations by the USP for microbial limit testing (3).

Tryptone and yeast extract provide nitrogenous and carbonaceous compounds, vitamin B complex and other growth nutrients. Dipotassium hydrogen phosphate provides buffering to the medium. During the first 24 hours, contaminating organisms are almost inhibited by tellurite, lithium chloride and high glycine content. The effect of inhibitors on *S.aureus* is reduced because of the presence of mannitol and glycine. Coagulase-positive staphylococci reduce potassium tellurite to metallic free tellurium and thus produce black colonies surrounded by yellow zones. This yellow colour is due to phenol red indicator that turns yellow in acidic condition due mannitol fermentation. If mannitol is not fermented, yellow zones are not formed. Also the colour of the medium around the colonies may even be a deeper red than normal due to utilization of the peptones in the medium. Prolonged incubation may result in growth of black coagulase-negative colonies.

Type of specimen

Clinical samples - pus samples, urine, wound samples; Food and dairy samples.

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (7,8).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (9).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further biochemical testing is required on colonies of pure culture for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to pink homogeneous free flowing powder

Gelling

Firm, comparable with 1.6% Agar gel.

Colour and Clarity of prepared medium

Red coloured clear to slightly opalescent gel forms in Petri plates.

Reaction

Reaction of 6.1% w/v aqueous solution at 25°C. pH : 7.2±0.2

pH

7.00-7.40

Cultural Response

Cultural characteristics observed with added 1% Potassium Tellurite solution (FD052), after an incubation at 35-37°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony	Mannitol fermentation
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 ⁴	inhibited	0%	-	-
<i>Proteus mirabilis</i> ATCC 25933	50-100	poor	10-20%	black	negative
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	luxuriant	≥50%	black with yellow halo	positive
<i>Staphylococcus epidermidis</i> ATCC 12228 (00036*)	50-100	Fair-good	30-40%	translucent to blackish	negative
<i>Escherichia coli</i> ATCC 8739 (00012*)	≥10 ⁴	inhibited	0%	-	-
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 6538 (00032*)	50-100	luxuriant	≥50%	black with yellow halo black with yellow halo	positive positive

Key : (*) Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (7,8).

Reference

1. FDA Bacteriological Analytical Manual, 2016, AOAC, Washington, D.C.
2. Vogel R. A. and Johnson M. J., 1960, Public Health Lab. 18:131.
3. United States Pharmacopeia, 2019. United States Pharmacopeial Convention, Inc., Rockville, Md.
4. Zebrovitz E., Evans J. B. and Niven C. F., 1955, J. Bacteriol., 70:686.
5. MacFaddin J. F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams & Wilkins, Baltimore, Md.
6. Curry A. S., Graf J. G. and McEwen G. M., (Eds.), 1993, CTFA Microbiology Guidelines, The Cosmetic, Toiletry and Fragrance Association, Washington, D.C.
7. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
8. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
9. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.


Revision : 04/2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

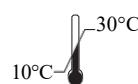

Storage temperature

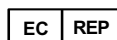

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.