

Haemophilus Test Agar Plate

Intended use

Recommended for the susceptibility testing of *Haemophilus influenzae* .

Composition**

Ingredients	g / L
HM infusion B from 300g #	2.000
Acicase ##	17.500
Yeast extract	5.000
Starch	1.500
Agar	17.000
Haematin Growth Supplement (FD117)	1 vial
NAD (2x7.50 mg)	15.000 mg
Haematin (2x7.50 mg)	15.000mg
Final pH (at 25°C)	7.4±0.2

**Formula adjusted, standardized to suit performance parameters

- Equivalent to Beef infusion from, ##-Equivalent to casein acid hydrolysate

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 cfu) aseptically on the plate.

Principle And Interpretation

Haemophilus species are nutritionally fastidious in nature. They require either exogenous hemin (X-factor) or NAD (V-factor) or both (1). Due to this reason, Mueller Hinton Agar, which is used for antimicrobial susceptibility of bacteria (2,3,4), cant be used for the antimicrobial susceptibility testing of *Haemophilus*. Haemophilus Test Agar Base, studied by Jorgensen et al (5,6) is used for the susceptibility testing of *Haemophilus influenzae*. This medium has similar composition as Mueller Hinton Agar, with the addition of yeast extract and added growth supplements. Haemophilus Test Agar Base is simple, transparent and posses minimum risk of antagonism of antimicrobial agents (5). Haemophilus Test Agar Base is also recommended by the United States National Committee for Clinical Laboratory Standards (NCCLS) for both dilution and disc diffusion assays (7). This medium scores over Mueller Hinton Agar with heamoglobin over clarity, thereby enabling proper visualization of inhibition zones. It also has low levels of the nucleotide thymidine, which allows testing of trimethoprim / sulphamethoxazole.

Haemophilus Test Agar Base contains HM infusion B from and acicase, which provide essential nutrients to the organisms. Yeast extract serves as a source of B complex vitamins. Starch acts as a protective colloid against toxic substances present in the medium. The surface of a Haemophilus Test Agar Base with added nutrients is inoculated either by using swab or by spreading the suspension. Antimicrobial discs i.e. paper discs impregnated with specific amount of antibiotics or other antimicrobial agents are placed on the surface of medium spaced properly. The plates are incubated in a CO₂ incubator and subsequently the inhibition zones around each disc are read. Comparing the zones of inhibition with the NCCLS standards, the determination as to whether the organism is susceptible, resistant or intermediate in its response to the antimicrobial substances is made (7).

Type of specimen

Isolated Microorganism from clinical samples.

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. The plates are to be incubated at 5-7% carbon dioxide at 35-37°C for 18-24 hours for appropriate results.
2. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
3. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile Haemophilus Test Agar in 90mm disposable plate with smooth surface and absence of black particles/cracks/bubbles

Colour of medium

Light amber coloured medium

Quantity of medium

25 ml of medium in 90 mm plate

pH

6.80 - 7.20

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed in 5-7% carbon dioxide after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Haemophilus influenzae</i> ATCC 49766	50-100	luxuriant	≥70%
<i>Enterococcus faecalis</i> ATCC 29212 (00187*)	50-100	good-luxuriant	≥70%
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant	≥70%
<i>Neisseria meningitidis</i> ATCC 13090	50-100	good-luxuriant	≥70%
<i>Staphylococcus aureus</i> subsp.aureus ATCC 25923 (00034*)	50-100	good-luxuriant	≥70%

(*) - Corresponding WDCM numbers

Storage and Shelf Life

Store between 2-8°C. Use before expiry period on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,5).

Reference

1. Barry A. L., Garcia F., and Thrupp L. D., 1970, Am. J. Clin. Pathol., 53 :149.
2. Bauer A. W., Kirby W. M., Sherris J. C. and Turck M., 1966, Am. J. Clin. Pathol. 45:493.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
4. Jorgensen J. H., Redding J. S., Maher L. A. and Howell A. W., 1987, J. Clin. Microbiol., 25:2105.
5. Jorgensen J. H., Howell A. W., and Maher L. A., J. Clin. Microbiol, 28:985.
6. Murray P. R., Baron J. H., Pfaller M. A., Jorgensen J. H. and Tenover F. C., (Ed.), 2003, Manual of Clinical Microbiology, 8th Ed., American Society for Microbiology, Washington, D.C.
7. Ryan K. J., Schoenkecht F. D., and Kirby W. M., 1970, Hospital Practice, 5:91.


Revision : 03/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



In vitro diagnostic
medical device

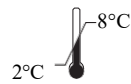

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



Do not use if
package is damaged

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.