

HiCrome™ UTI Agar Plate

MP1353

Intended use

Recommended for presumptive identification and confirmation of microorganisms mainly causing urinary tract infections, can also be used for testing water, food, environmental and other clinical samples.

Composition**

Ingredients	g / L
Peptone, special	15.000
Chromogenic mixture	2.450
Agar	15.000
Final pH (at 25°C)	6.8±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

Principle And Interpretation

Urinary tract infections are bacterial infections affecting parts of urinary tract. The common symptoms of urinary tract infection are urgency and frequency of micturition, with associated discomfort or pain. The common condition is cystitis, due to infection of the bladder with a uropathogenic bacterium, which most frequently is *Escherichia coli*, but sometimes *Staphylococcus saprophyticus* or especially in hospital-acquired infections, *Klebsiella* species, *Proteus mirabilis*, other coliforms, *Pseudomonas aeruginosa* or *Enterococcus faecalis* (1). HiCrome™ UTI Agar is formulated on basis of work carried out by Pezzlo (2) Wilkie et al (3), Friedman et al (4), Murray et al (2), Soriano and Ponte (5) and Merlino et al (6). These media are recommended for the detection of urinary tract pathogens where HiCrome™ UTI Agar has broader application as a general nutrient agar for isolation of various microorganisms. It facilitates and expedites the identification of some gram-negative bacteria and some gram-positive bacteria on the basis of different contrasted colony colours produced by reactions of genus or species specific enzymes with two chromogenic substrates. The chromogenic substrates are specifically cleaved by enzymes produced by *Enterococcus* species, *E.coli* and coliforms. Presence of amino acids like phenylalanine and tryptophan from peptones helps for detection of tryptophan deaminase activity, indicating the presence of *Proteus* species, *Morganella* species and *Providencia* species.

One of the chromogenic substrate is cleaved by β -glucosidase possessed by Enterococci resulting in formation of blue colonies. *E.coli* produce pink colonies due to the enzyme β -D-galactosidase that cleaves the other chromogenic substrate. Further confirmation of *E.coli* can be done by performing the indole test. Coliforms produce purple coloured colonies due to cleavage of both the chromogenic substrate. Colonies of *Proteus*, *Morganella* and *Providencia* species appear brown because of tryptophan deaminase activity. Peptone special provides nitrogenous, carbonaceous compounds, long chain amino acids, vitamins and other essential growth nutrients. This medium can be made selective by supplementation with antibiotics for detecting microorganisms associated with hospital borne infections.

Type of specimen

Clinical samples : urine, faeces, etc.; Food samples , Water samples.

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,7).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (3). For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (8).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Since it is an enzyme-substrate based reaction, the intensity of colour may vary with isolates.
2. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium
3. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile HiCrome™ UTI Agar in 90mm disposable plates with smooth surface and absence of black particles/cracks/bubbles

Colour of medium

Light amber coloured

Quantity of medium

25ml of medium in 90mm disposable plate

pH

6.60-7.00

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 16-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of Colony
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	>=70%	Purple to magenta
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	luxuriant	>=70%	blue-green (small)
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	50-100	luxuriant	>=70%	blue to purple, mucoid
<i>Proteus mirabilis</i> ATCC 12453	50-100	luxuriant	>=70%	light brown
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	luxuriant	>=70%	colourless (greenish pigment may be observed)
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	luxuriant	>=70%	cream to golden yellow

Key : (*) Corresponding WDCM numbers.

Storage and Shelf Life

On receipt store between 2-8°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,7).

Reference

1. Collee J. G., Fraser A. G., Marmion B. P., Simmons A., (Eds.), Mackie and McCartney, Practical Medical Microbiology, 1996, 14th Edition, Churchill Livingstone.
2. Murray P., Traynor P. Hopson D., 1992, J. Clin. Microbiol., 30:1600- 1601.
3. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.

- 4.Friedman M. P. et al, 1991, J. Clin. Microbiol., 29:2385-2389.
- 5.Soriano F., Ponte C., 1992, J. Clin. Microbiol., 30:3033-3034.
- 6.Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
- 7.Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
- 8.Lipps WC, Braun-Howland EB, Baxter TE,eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.


Revision : 03/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



IVD *In vitro* diagnostic
medical device

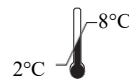

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.