



# Technical Data

## MacConkey Agar w/ CV, NaCl and 0.15% Bile salts

M081A

### Intended Use:

Recommended to identify *Enterobacteriaceae* in the presence of coliforms and lactose nonfermenters from water, sewage, food products etc.

### Composition\*\*

Ingredients	Gms / Litre
Peptone	20.000
Lactose	10.000
Sodium chloride	5.000
Bile salts	1.500
Neutral red	0.050
Crystal violet	0.001
Agar	15.000
Final pH ( at 25°C)	7.2±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 51.55 grams of medium in 1000 ml purified / distilled water. Heat to boiling with gentle swirling to dissolve the agar completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Avoid overheating. Cool to 45 - 50°C. Mix well and pour into sterile Petri plates. The surface of the medium should be dry when inoculated.

### Principle And Interpretation

MacConkey Agar is the earliest selective and differential medium for cultivation of enteric microorganisms from a variety of clinical specimens (5,6). Subsequently MacConkey Agar and Broth have been recommended for use in microbiological examination of foodstuffs (7) and for direct plating / inoculation of water samples for coliform counts (2). These media are also accepted by the Standard Methods for the Examination of Milk and Dairy Products (8). Original medium contains protein, bile salts, sodium chloride and two dyes.

The selective action of this medium is attributed to crystal violet and bile salts, which are inhibitory to most species of gram-positive bacteria. Gram-negative bacteria usually grow well on the medium and are differentiated by their ability to ferment lactose. Lactose fermenting strains grow as red or pink and may be surrounded by a zone of acid precipitated bile. The red colour is due to production of acid from lactose, absorption of neutral red and a subsequent colour change of the dye when the pH of medium falls below 6.8. Lactose non-fermenting strains, such as *Shigella* and *Salmonella* are colourless and transparent and typically do not alter appearance of the medium. *Yersinia enterocolitica* may appear as small, non-lactose fermenting colonies after incubation at room temperature.

Peptone serves as a source of carbon, nitrogen, long chain amino acids, vitamins and other essential growth nutrients. Lactose serves as the carbon source by being the fermentable carbohydrate. Bile salts and crystal violet serves to make the medium selective by inhibiting accompanying gram-positive bacteria. Neutral red is the pH indicator dye while sodium chloride maintains the osmotic equilibrium of the medium.

### Type of specimen

Clinical samples: faeces, urine and other pathological material; Food and dairy samples; water samples.

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (1,7,8).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (2).

After use, contaminated materials must be sterilized by autoclaving before discarding.

## Warning and Precautions

In Vitro diagnostic use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

## Limitations

1. Though the medium is recommended for selective isolation, further biochemical and serological testing must be carried out for complete identification.
2. The surface of the medium should be dry when inoculated.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Light yellow to pink homogeneous free flowing powder

### Gelling

Firm comparable with 1.35% Agar gel.

### Colour and Clarity of prepared medium

Red with purplish tinge coloured clear to slightly opalescent gel forms in Petri plates.

pH of 4.95% w/v aqueous solution at 25°C (after sterilization).

### pH

7.00-7.40

### Cultural Response

Growth Promotion is carried out in accordance with the harmonized method of ICH (USP/EP/BP/JP). Cultural response was observed after an incubation at 30-35°C for 18-72 hours. Recovery rate is considered as 100% for bacteria growth on Soybean Casein Digest Agar.

### Growth promoting properties

Growth of microorganism comparable to that previously obtained with previously tested and approved lot of medium occurs at the specified temperature for not more than the shortest period of time specified inoculating 100 cfu (at 30-35°C for 18 hours).

### Indicative properties

Colonies are comparable in appearance and indication reaction to those previously obtained with previously tested and approved lot of medium occurs for the specified temperature for a period of time within the range specified inoculating  $\leq 100$  cfu (at 30-35°C for 18-72 hours).

Organism	Inoculum (CFU)	Growth	Observed Lot value (CFU)	Recovery	Colour of colony	Incubation temperature
<b>Growth Promoting + Indicative</b>						
<i>Escherichia coli</i> ATCC 8739 (00012*)	50 -100	luxuriant	25 -100	$\geq 50$ %	pink-red with bile precipitate	18 -72 hrs
<b>Additional Microbiological testing</b>						
<i>Escherichia coli</i> ATCC 25922 (00013*)	50 -100	luxuriant	25 -100	$\geq 50$ %	pink to red with bile precipitate	18 -24 hrs
<i>Escherichia coli</i> NCTC 9002	50 -100	luxuriant	25 -100	$\geq 50$ %	pink to red with bile precipitate	18 -24 hrs
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50 -100	luxuriant	25 -100	$\geq 50$ %	pink to red	18 -24 hrs
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50 -100	none-poor	0 -10	0 -10 %	colourless to pale pink	18 -24 hrs
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs

<i>Staphylococcus aureus</i> <i>subsp. aureus</i> ATCC 6538 (00032*)	$\geq 10^4$	inhibited	0	0 %		$\geq 24$ hrs
<i>Staphylococcus aureus</i> <i>subsp. aureus</i> ATCC 25923 (00034*)	$\geq 10^4$	inhibited	0	0 %		$\geq 24$ hrs
<i>Salmonella</i> Enteritidis ATCC 50 -100 13076 (00030*)		luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Salmonella</i> Paratyphi A ATCC 9150	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Salmonella</i> Paratyphi B ATCC 8759	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Salmonella</i> Typhi ATCC 6539	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Salmonella</i> Abony NCTC 6017 (00029*)	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Proteus vulgaris</i> ATCC 13315	50 -100	luxuriant	25 -100	$\geq 50$ %	colourless	18 -24 hrs
<i>Shigella flexneri</i> ATCC 12022 (00126*)	50 -100	fair to good	15 -40	30 -40 %	colourless	18 -24 hrs
<i>Staphylococcus epidermidis</i> ATCC 12228 (00036*)	$\geq 10^4$	inhibited	0	0 %		$\geq 24$ hrs
<i>Corynebacterium</i> <i>diphtheriae</i> type <i>gravis</i>	$\geq 10^4$	inhibited	0	0 %		$\geq 24$ hrs

Key : (#) *Enterobacter aerogenes* ATCC 13048 (\*) Corresponding WDCM numbers

## Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,4).

## Reference


1. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
2. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2<sup>nd</sup> Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
5. MacConkey, 1900, The Lancet, ii:20.
6. MacConkey, 1905, J. Hyg., 5:333.
7. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
8. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C. 7.



In vitro diagnostic medical device



CE Marking

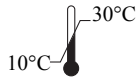

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,  
23 Vadhani Industrial Estate,  
LBS Marg, Mumbai-86, MS, India



CE Partner 4U ,Esdoornlaan 13, 3951  
DB Maarn The Netherlands,  
[www.cepartner4u.eu](http://www.cepartner4u.eu)

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.