



Nutrient Broth

M002

Intended use

Nutrient Broth is used for the general cultivation of less fastidious microorganisms, can be enriched with blood or other biological fluids.

Composition**

Ingredients	Gms / Litre
Peptone	5.000
Sodium chloride	5.000
HM peptone B [#]	1.500
Yeast extract	1.500
Final pH (at 25°C)	7.4±0.2

**Formula adjusted, standardized to suit performance parameters

- Equivalent to Beef extract

Directions

Suspend 13.0 grams in 1000 ml distilled water. Heat, if necessary, to dissolve the medium completely. Dispense into tubes or flasks as desired. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle And Interpretation

Nutrient media are basic culture media used for maintaining microorganisms, cultivating fastidious organisms by enriching with serum or blood and are also used for purity checking prior to biochemical or serological testing (6,7). Nutrient Broth has the formula originally designed for use in the Standard Method for Examination of Water and Waste water. It is one of the several non-selective media useful in routine cultivation of microorganisms (1, 2). It can be used for the cultivation and enumeration of bacteria which are not particularly fastidious. Addition of different biological fluids such as horse or sheep blood, serum, egg yolk etc. makes it suitable for the cultivation of related fastidious organisms.

Peptone, HM peptone B and yeast extract provide the necessary nitrogen compounds, carbon, vitamins and also some trace ingredients necessary for the growth of bacteria. Sodium chloride maintains the osmotic equilibrium of the medium.

Type of specimen

Clinical samples - Blood ; Food and dairy samples ; Water samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,7).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (1,2,8).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards.(3)

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

Limitations :

This medium is general purpose medium and may not support the growth of fastidious organisms.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Colour and Clarity of prepared medium

Light yellow coloured clear to slightly opalescent solution

Reaction

Reaction of 1.3% w/v aqueous solution at 25°C. pH : 7.4±0.2

pH

7.20-7.60

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-48 hours.

Cultural Response

Organism	Inoculum (CFU)	Growth
Cultural Response		
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	good-luxuriant
<i>Salmonella Typhi</i> ATCC 6539	50-100	good-luxuriant
<i>Staphylococcus aureus</i> ATCC 25923 (00034*)	50-100	good-luxuriant
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store below 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label.

Product performance is best if used within stated expiry period. .

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference

1. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
2. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2001, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
3. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
4. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.

6. Lapage S., Shelton J. and Mitchell T., 1970, 'Methods in Microbiology', Norris J. and Ribbons D., (Eds.), Vol. 3A, Academic Press, London.
7. MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Lippincott, Williams and Wilkins, Baltimore.
8. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 03 / 2018



In vitro diagnostic medical device



CE Marking

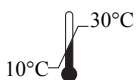

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.