



Hoyle Medium Base

M015

Intended use

Hoyle Medium Base is a highly selective medium used for the isolation and differentiation of *Corynebacterium diphtheriae* types.

Composition**

Ingredients	Gms / Litre
Peptone	10.000
HM extract B#	10.000
Sodium chloride	5.000
Agar	15.000
Final pH (at 25°C)	7.8±0.2

**Formula adjusted, standardized to suit performance parameters

- Equivalent to Beef extract

Directions

Suspend 40.0 grams in 940 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C and aseptically add 50 ml of laked blood and 10 ml of 3.5% Potassium Tellurite Solution (FD047). Mix well and pour into sterile Petri plates.

Principle And Interpretation

The most common disease caused by *Corynebacterium diphtheriae* is diphtheria, an acute communicable disease manifested by both local infection of the upper respiratory tract and the systemic effects of the toxin, which are most notable in the heart and peripheral nerves (1). Hoyle Medium Base, formulated by Hoyle (2), is the modification of the original formulation of Neill, for the isolation and differentiation of *C.diphtheriae*. This medium is not inhibitory to some mitis types of *Corynebacterium*, as the original formulation.

Peptone and HM extract supply essential growth nutrients. Potassium tellurite is a selective agent, which inhibits most of the normal flora of the upper respiratory tract except *Corynebacterium*. Hoyle's Medium is a highly selective medium and should be used in conjunction with a non-selective media such as Loeffler Serum Medium (M537) and Blood Agar Base (M089) with 10% horse blood (3).

C.diphtheriae are usually present in small numbers permitting the formation of well isolated colonies. So, inoculation is done by directly rubbing the swab over the entire surface of the medium. Incubation should be carried out till 72 hours if the results are negative. To study the morphology, gentian violet staining is done. To demonstrate the characteristic morphology and staining reactions of *C. diphtheriae* by Neissers or Alberts method, it is advisable to use colonies from Loeffler Medium. The toxigenicity of *C.diphtheriae* strains can be determined by Eleks (4) method.

Type of specimen

Clinical samples - Blood ; Food and dairy samples ; Water samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

Limitations :

This medium is general purpose medium and may not support the growth of fastidious organisms.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Basal Medium: Amber coloured, clear to slightly opalescent gel . After Addition of blood & Tellurite : Brownish red coloured opaque gel forms in Petri plates

Reaction

Reaction of 4.0% w/v aqueous solution at 25°C. pH : 7.8±0.2

pH

7.60-8.00

Cultural Response

M015: Cultural characteristics observed with added 50 ml of laked blood and tellurite solution, after an incubation at 35-37°C for 18-24 hours .

Organism	Inoculum (CFU)	Growth	Recovery	Colony characteristics
Cultural Response <i>Bacillus subtilis</i> ATCC 6633 (00003*)	>=10 ³	inhibited	0%	
<i>C. diphtheriae</i> type <i>intermedius</i> 14779	50-100	good-luxuriant	>=50%	grey colonies with darker centers
<i>Corynebacterium diphtheriae</i> type <i>mitis</i>	50-100	good-luxuriant	>=50%	grey colonies with shining surface
<i>Escherichia coli</i> ATCC 25922 (00013*)	>=10 ³	inhibited	0%	
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	good-luxuriant	>=50%	black minute colonies

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label.

Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference


- 1.Murray P. R., Baron J. H., Pfaller M. A., Jorgensen J. H. and Tenover F. C., 8th Ed., American Society for Microbiology, Washington, D.C. (Ed.), 2003, Manual of Clinical Microbiology,
- 2.Hoyle I., 1941, Lancet., 1:175.
- 3.MacFaddin J. F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams and Wilkins, Baltimore.
- 4.Elek S. D., 1948, Brit. Med. A1:493.
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.



In vitro diagnostic medical device



CE Marking

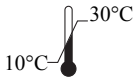

Storage temperature




Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
B /4-6 , MIDC, Dindori, Nashik MH

www.himedialabs.com

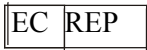

CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.