



Maximum Recovery Diluent

M1030

Intended use

Maximum Recovery Diluent is a protective and isotonic medium used for maximal recovery of microorganisms from a variety of sources.

Composition**

Ingredients	Gms / Litre
Peptone	1.000
Sodium chloride	8.500
Final pH (at 25°C)	7.0±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 9.5 grams in 1000 ml purified/distilled water. Mix well and dispense into tubes or flasks as desired. Heat if necessary to dissolve the medium completely. Cool to 45-50°C. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle And Interpretation

Maximum Recovery Diluent is formulated as recommended by ISO Committee (1) for use as an isotonic diluent. Standard methods for the examination of foodstuffs require sample dilution to be carried out accurately to estimate the number of microorganisms.

Maximum Recovery Diluent combines protective effect of peptone (6) with the osmotic balance of physiological saline (4). The low concentration of peptone helps to maintain the organisms for 1-2 hours of dilution without multiplication. The isotonic property of the diluent ensures the recovery of organisms from various sources, which may be vulnerable in distilled water or aqueous suspensions.

Put 10 gm of test sample into a sterile blender jar and add 90 ml of sterile Maximum Recovery Diluent. Operate the blender at 15,000 to 20,000 revolutions per minute. Transfer 1 ml of it to 9 ml of sterile diluent within 15 minutes and mix well. It will be 10⁻¹ dilution. Appropriate serial dilutions can be prepared using same diluent and counts obtained by spread plate or pour plate technique. Use a positive test sample divided between new and previous diluent. Carry out duplicate tests as described in technique and look for equivalent yields of organisms between the diluent batches. Incubate the tubes with test organisms. At time zero and after 30 minutes at room temperature, subculture a loopful (0.01 ml) onto Soyabean Casein Digest Agar (M290) with 5% v/v sheep blood using streak technique. Incubate plates at 35 ± 2°C for 18-24 hours.

Type of specimen

Clinical samples- urine samples , stool samples, Food samples.

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (2,3).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Due to nutritional variations some organisms may show poor growth.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

White to pale yellow homogeneous free flowing powder

Colour and Clarity of prepared medium

Light yellow coloured clear solution without any precipitate

Reaction

Reaction of 0.95% w/v aqueous solution at 25°C. pH : 7.0±0.2

pH

6.80-7.20

Cultural Response

Cultural characteristics observed on Soyabean Casein Digest Agar (M290), after an incubation at 35-37°C for 18-24 hours of cultures suspended in Maximum Recovery Diluent for 30 minutes.

Organism	Inoculum (CFU)	Recovery (after 30 minutes)
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	no change in numbers
<i>Staphylococcus aureus</i> subsp <i>aureus</i> ATCC 25923 (00034*)	50-100	no change in numbers

Key : (*) Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 15-25°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label.

Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (2,3).

Reference

1. International Organization for the Standardization (ISO), ISO/DIS 6649.
2. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
3. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
4. Patterson J. W. and Cassells J. A., 1963, J. Appl. Bacteriol., 26:493.
5. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Straker R. P. and Stokes J. L., 1957, Appl. Microbiol., 5:21.


Revision : 03 / 2019



In vitro diagnostic medical device



CE Marking

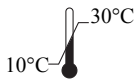

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U, Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.