



## Hemorrhagic Coli (HC) Agar

M1158

### Intended use

Hemorrhagic coli (HC) Agar is recommended for isolation and enumeration of *Escherichia coli* with an enzyme labelled monoclonal antibody.

### Composition\*\*

Ingredients	Gms / Litre
Tryptone	20.000
Sorbitol	20.000
Sodium chloride	5.000
Bile salts	1.120
Bromocresol purple	0.015
Agar	15.000
Final pH ( at 25°C)	7.2±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 61.13 grams in 1000 ml purified/ distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

Hemorrhagic colitis is a type of gastroenteritis in which certain strains of the bacterium *Escherichia coli* infect the large intestine and produce a toxin that causes bloody diarrhea and other serious complications. *Escherichia coli* O157:H7 was recognized as the cause of hemorrhagic colitis (3). Outbreaks can be caused by eating undercooked beef, especially ground beef, or by drinking unpasteurized milk or juice. HC Agar is used for isolation and enumeration of *E.coli*, the primary strain causing hemorrhagic colitis (4,5).

Tryptone in the medium is a source of carbon, nitrogen, vitamins and minerals. Sodium chloride maintains the osmotic balance of the medium. *E.coli* O157:H7 does not ferment sorbitol and produces colourless colonies. Bromocresol purple acts as a pH indicator. Bile salts inhibit accompanying gram-positive bacteria.

Homogenize 10 grams of test sample in 90 ml Peptone Water (M028) (prepare 1:100 dilutions if counts are expected to be high). Pipette 1 ml aliquots through disposable 100 µm pre filter and add to 10 ml Peptone Water (M028) filtered through Hydrophobic Grid Membrane Filters (HGFM) (3). Lay the filters onto Hemorrhagic Coli (HC) Agar and incubate at 43°C for 16-20 hours. Replicate colony growth onto other HGMFs using HGFM replicator. Incubate replicates on HC Agar at 43°C for 16-20 hours and test original filters with conjugated antibody (4).

### Type of specimen

Clinical samples: Blood

### Specimen Collection and Handling

For clinical samples follow appropriate techniques or handling specimens as per established guidelines (1,2).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations :

1. Some strains may show poor growth due to nutritional variation.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Cream to yellow homogeneous free flowing powder

### Gelling

Firm, comparable with 1.5% Agar gel

### Colour and Clarity of prepared medium

Purple coloured, clear to slightly opalescent gel forms in Petri plates

### Reaction

Reaction of 6.11% w/v aqueous solution at 25°C. pH : 7.2±0.2

### pH

7.00-7.40

### Cultural Response

Cultural characteristics observed after an incubation at 43°C for 16-20 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Escherichia coli</i> O157:H7	50-100	good-luxuriant	≥50%
<i>Proteus mirabilis</i> ATCC 25933	50-100	good	40-50%

## Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (1,2).

## Reference

1. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
2. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
3. Riley L. W., Remis R. S., Helgerson S. D., 1983, N. Engl. J. Med., 308:681.
4. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
5. Todd et al, 1988, Appl. Environ. Microbiol., 54:2536.


Revision : 04 / 2019



In vitro diagnostic medical device



CE Marking

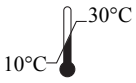

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,  
23 Vadhani Industrial Estate,  
LBS Marg, Mumbai-86, MS, India



CE Partner 4U .Esdoornlaan 13, 3951  
DB Maarn The Netherlands,  
[www.cepartner4u.eu](http://www.cepartner4u.eu)

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.