



HiCrome™ L.mono Rapid Differential Agar Base

M1924

Intended Use:

Recommended for the rapid identification and differentiation of *Listeria monocytogenes* from other *Listeria* species based on rhamnose fermentation and PIPLC activity.

Composition**

Ingredients	Gms / Litre
Peptone special	23.000
Tryptone	10.000
Soya peptone	2.000
Sodium chloride	4.000
Lithium chloride	5.000
Chromogenic mixture	1.160
Rhamnose	10.000
Phenol red	0.120
Agar	15.000
Final pH (at 25°C)	7.4±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 35.14 grams in 470 ml purified /distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Aseptically add sterile contents of 1 vial of L. mono Enrichment Supplement I (FD214) and sterile rehydrated contents of 1 vial of HiCrome™ Listeria Selective Supplement (FD181). Mix well and pour into sterile Petri plates.

Principle And Interpretation

Listeria monocytogenes is a gram-positive foodborne human pathogen responsible for serious infections in pregnant women that may ultimately result in abortion, stillbirth, birth of a child with neonatal listeriosis and meningitis or primary bacteremia in adults and juveniles. The pathogenicity of *Listeria ivanovii* for humans is uncertain (6). Since *L. monocytogenes* and *L.*

innocua have similar biochemical properties, they cannot be differentiated on traditional media (PALCAM, Oxford). This medium is based on the specific chromogenic detection of β-glucosidase activity, rhamnose fermentation and PIPLC activity. *Listeria* species hydrolyse the purified chromogenic substrate in the medium giving blue coloured colonies. Since β-glucosidase activity is specific for *Listeria* species, other organisms cannot utilize the chromogenic substrate and therefore give white colonies. Differentiation between *Listeria* species is based on the property of rhamnose fermentation and PIPLC

activity. The colonies of *L. monocytogenes* appear bluish green with a yellow halo (rhamnose positive) while the colonies of *L. ivanovii* appear bluish green without a yellow halo (Rhamnose negative) (3,4). The differentiation of *L. monocytogene* and *L. innocua* is based on PIPLC phosphatidylinositol-specific phospholipase C activity. Phospholipase C enzyme hydrolyses the purified substrate (FD214) added to the medium resulting in an opaque halo around *Listeria monocytogenes* colonies. *L. ivanovii* also demonstrates PIPLC activity however since it does not ferment rhamnose it can be easily distinguished from *L. monocytogenes* (5,6).

Peptone special, tryptone and soya peptone provide nitrogenous, carbonaceous substances, long chain amino acids, vitamin B complex and other essential growth nutrients. Rhamnose is the fermentable carbohydrate with phenol red as an indicator. Sodium chloride maintains the osmotic equilibrium. The added lithium chloride and HiCrome™ Listeria Selective Supplement (FD181) inhibit growth of most gram-positive bacteria, gram-negative bacteria, yeasts and moulds. Phospholipase C enzyme hydrolyses the purified substrate (FD214) added to the medium resulting in an opaque halo around *Listeria monocytogenes* colonies demonstrating PIPLC activity.

Type of specimen

Clinical sample; Food samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (1,2).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (7,8). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Due to variable nutritional requirements, some strains may show poor growth on this medium.
2. Slight colour variation may be observed depending upon the utilization of the substrate by the organisms.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to pink homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel.

Colour and Clarity of prepared medium

Red coloured, opalescent gel forms in Petri plates

Reaction

Reaction of 7.03% w/v aqueous solution at 25°C. pH : 7.4±0.2

pH

7.20-7.60

Cultural Response

Cultural characteristics observed w/added HiCrome™ Listeria Selective Supplement (FD181) and L.mono Enrichment supplement I (FD214), after an incubation at 35-37°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony	Rhamnose fermentation	PIPLC Activity
<i>Listeria innocua</i> ATCC 33090 (00017*)	50-100	luxuriant	≥50%	bluish green	positive reaction, (yellow background)	negative reaction
<i>Listeria ivanovii</i> subsp. <i>ivanovii</i> serovar 5 ATCC 19119 (00018*)	50-100	luxuriant	≥50%	bluish green	negative reaction	positive, opaque halo around the colony exhibiting phosphatidyl inositol specific phospholipase activity.
<i>Listeria monocytogenes</i> ATCC 19118	50-100	luxuriant	≥50%	bluish green	positive reaction, (yellow background)	positive, opaque halo around the colony exhibiting phosphatidyl inositol specific phospholipase activity.
<i>Listeria monocytogenes</i> ATCC 13932 (00021*)	50-100	luxuriant	≥50%	bluish green	positive reaction, (yellow background)	positive, opaque halo around the colony exhibiting phosphatidyl inositol specific phospholipase activity.
<i>Listeria monocytogenes</i> ATCC 35152 (00109*)	50-100	luxuriant	≥50%	bluish green	positive reaction, (yellow background)	positive, opaque halo around the colony exhibiting phosphatidyl inositol specific phospholipase activity.

<i>Bacillus subtilis subsp. spizizenii</i> ATCC 6633 (00003*)	$\geq 10^4$	Inhibited	0%
<i>Escherichia coli</i> ATCC 25922 (00013*)	$\geq 10^4$	inhibited	0%
<i>Escherichia coli</i> ATCC 8739 (00012*)	$\geq 10^4$	inhibited	0%
<i>Enterococcus faecalis</i> ATCC 19433 (00009*)	$\geq 10^4$	inhibited	0%
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	$\geq 10^4$	inhibited	0%

Key : (*) Corresponding WDCM numbers.

Storage and Shelf Life

Store between 15-25°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (1,2).

Reference

1. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
2. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
3. Mengaud J., Braun-Breton C. and Cossart P., (1991), Molecular Microbiology, 5(2): 367-372.P
4. Notermans S.H. and Dufrenne J., (1991), Applied and Environmental Microbiology, 57(09): 2666-70.
5. Ottaviani F., Ottaviani M., and Agosti M. (1997 a), Industrie Alimentari 36, 1-3.
6. Ottaviani F., Ottaviani M., and Agosti M. (1997 b), Quimper Froid Symposium Proceedings p. 6, A.D.R.I.A. Quimper,France, 16-18 June 1997.
7. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
8. Schlech WF, Lavigne PM, Bortolussi RA, et al.(January 1983). "Epidemic listeriosis-evidence for transmission by food". N. Engl. J.Med. 308(4): 203–6. doi:10.1056.


Revision : 03/2020



In vitro diagnostic medical device



CE Marking

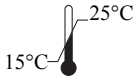

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.