



## HM Peptone B Agar

M806

### Intended use

HM Peptone B Agar is used as a general purpose nutrient medium which can support growth of not particularly fastidious bacteria

### Composition\*\*

Ingredients	Gms / Litre
Peptone	10.000
HM peptone B#	3.000
Sodium chloride	5.000
Agar	15.000
Final pH ( at 25°C)	7.6±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# Equivalent to Beef extract

### Directions

Suspend 33.0 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour the medium in sterile Petri plates.

### Principle And Interpretation

The majority of organisms to be studied in medical bacteriology are either pathogens or commensals of the human body, and in order to obtain suitable growth the artificial culture medium should provide nutrients and a pH (about 7.2) approximating to those of the tissues and body fluids. For routine purposes many of these nutrients are supplied by aqueous extracts of meat and peptone, which is a product of the digestion of protein (1).

Meat Extract Agar B can be used as a general-purpose nutrient medium and is also recommended for preparation of pure culture of *Candida* species for carrying out fermentation studies (2). Meat Extract B Agar is a non-selective nutrient medium containing meat extract B and peptone as a source of nitrogen and carbon and sodium chloride as a source of electrolytes.

### Type of specimen

Clinical samples - Blood ; Food and dairy samples ; Water samples

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (1,2,8).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards.(3)

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

### Limitations :

This medium is general purpose medium and may not support the growth of fastidious organisms.

## Quality Control

### Appearance

Cream to yellow homogeneous free flowing powder

### Gelling

Firm, comparable with 1.5% Agar gel

### Colour and Clarity of prepared medium

Yellow coloured, clear to slightly opalescent gel forms in Petri plates

### Reaction

Reaction of 3.3% w/v aqueous solution at 25°C. pH : 7.6±0.2

### pH

7.40-7.80

### Cultural Response

M806: Cultural characteristics observed after an incubation at 35-37°C for 18-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	luxuriant	≥70%
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	≥70%
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	luxuriant	≥70%
<i>Salmonella Typhi</i> ATCC 6539	50-100	luxuriant	≥70%
<i>Staphylococcus aureus</i> ATCC 25923 (00034*)	50-100	luxuriant	≥70%

Key : \*Corresponding WDCM numbers.

## Storage and Shelf Life

Store below 30°C in tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label. Store below 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period. .

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,7).

## Reference


1. Collee J. G., Fraser A. G., Marimon B. P., Simmons A., (Eds.), 1996, Mackie and McCartney Practical Medical Microbiology, 14th Ed., Churchill Livingstone.
2. Finegold S. M. and Baron E. J., (Ed.), Bailey and Scott's Diagnostic Microbiology, 1986, 7th Edition, The C.V. Mosby Company, St. Louis.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.



In vitro diagnostic medical device



CE Marking

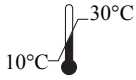

Storage temperature




Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,  
B /4-6 , MIDC, Dindori, Nashik MH

[www.himedialabs.com](http://www.himedialabs.com)

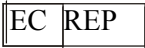

CE Partner 4U ,Esdoomlaan 13, 3951  
DB Maarn The Netherlands,

[www.cepartner4u.eu](http://www.cepartner4u.eu)

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.