

HIMEDIA®

For Life is Precious

*Unleash
the Power of
Automation*

Introducing HiPurA® Endotoxin-Free Plasmid DNA Purification Kits

HiGenoMB
Unzipping Genes

Insta NX® Mag24

Insta NX® Mag16



Single sample/
Pre-filled Cartridge



Pre-filled Plates



Single Pre-filled Cartridge
for processing Higher
sample volume



Achieve Seamless Extraction with HiPurA® Automated Extraction Kits



Miniprep (MB513PC16)

Yield: upto 100 ng/ μ l,
Purity ~1.8



Midiprep (MB514PC-24)

Yield: upto 2000 ng/ μ l,
Purity ~1.8



Maxiprep (MB515PC-24)

Yield: upto 4000 ng/ μ l,
Purity ~1.8

Endotoxin-free



Avoid Consumables
Wastage



Fast and Efficient



Convenient



Contamination
Free

Check out our
**Whitepapers on Endotoxin-Free
Plasmid DNA** for more information.



MB515PC-24



MB514PC-24

Order your Endotoxin-Free Plasmid DNA Purification Kit today!

Miniprep Insta NX® Mag16^{plus}


High quality
sequencing
grade Plasmid
DNA in 30 mins

Sample volume:
1mL (with
respect to O.D.)

Yield:
upto 100 ng/μl

Purity
(A_{260}/A_{280} ratio):
~ 1.8-1.9

MB513PC16 (Endotoxin free)

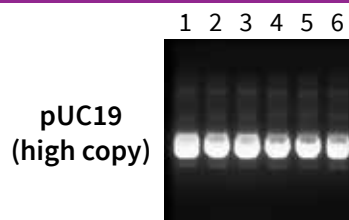

Lanes	Protocol	Conc. (ng/μl)	Purity (A_{260}/A_{280})
1	HiMedia	155.15	1.86
2		163.30	1.87
3		163.25	1.86
4	Competitor Q	81.40	1.92
5		77.35	1.92
6		79.55	1.91

Midiprep Insta NX® Mag24


High quality
sequencing
grade Plasmid
DNA in 60 mins

Sample volume:
150mL (with
respect to O.D.)

Yield:
1000-2000 ng/μl

Purity
(A_{260}/A_{280} ratio):
~ 1.8-1.9

MB514PC-24 (Endotoxin-free)

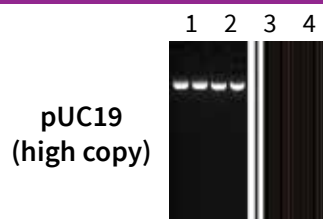

Lanes	Protocol	Conc. (ng/μl)	Purity (A_{260}/A_{280})
1	HiMedia	1916.50	1.91
2		2085.90	1.91
3	Competitor Q	815.15	1.93
4		899.25	1.93

Maxiprep Insta NX® Mag24


High quality
sequencing
grade Plasmid
DNA in 60 mins

Sample volume:
250mL (with
respect to O.D.)

Yield:
2500-4000 ng/μl

Purity
(A_{260}/A_{280} ratio):
~ 1.8-1.9

MB515PC-24 (Endotoxin free)

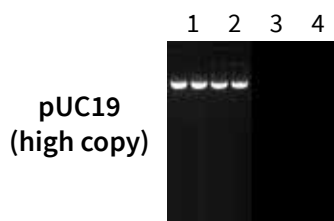

Lanes	Protocol	Conc. (ng/μl)	Purity (A_{260}/A_{280})
1	HiMedia	3291.10	1.91
2		3620.50	1.91
3	Competitor Z	1858.60	1.90
4		1737.00	1.92



For Life is Precious

HiMediaLaboratories™

www.himedialabs.com

CORPORATE OFFICE

Plot No. C40, Road No.21Y, MIDC, Wagle Industrial Estate, Thane (West) - 400604, Maharashtra, India.

Tel : +91-22-6147 1919 / 6116 9797 / 6903 4800 | Fax : +91-22-6147 1920

Email : mb@himedialabs.com



Literature Code : TL902_01/Plasmid DNA Extraction/0425