


Bacterial Endotoxin Assay Kit

Harmonized Gel Clot Limit Method for Testing of Bacterial Endotoxin
as per USP, EP, JP and ICH Guidelines

AccuBET™

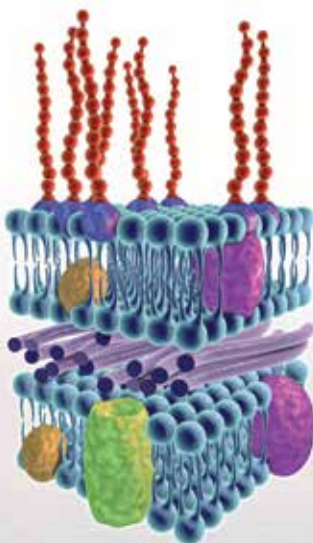

Popular in 15 Countries



No Lysate Reconstitution



Ampoule is the reaction Tube



Eliminate Reagent Wastage

New
Format

Single-Test / Ampoule Format

Gel - Clot Assay




MicroBiology™
Understanding Microbes

HIMEDIA®
For Life is Precious

Bacterial Endotoxin Assay Kit

Harmonized Gel-Clot Limit Method for Testing of Bacterial Endotoxin
as per USP, EP, JP and ICH Guidelines

ORIGIN OF ENDOTOXIN

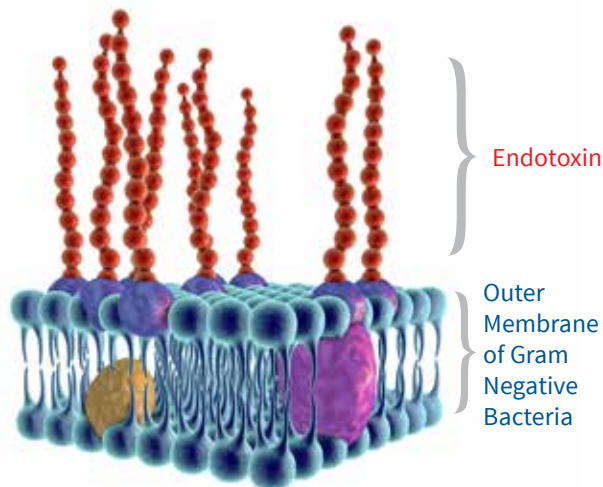

POTENTIAL USER SEGMENTS :

- Pharmaceuticals
- Vaccines
- Biologicals
- Hospitals
- Medical Device Manufacturers

ADVANTAGES :

- Single-test format available for all lysate sensitivities
- Ampoule packing reduces need for pyrogen free tubes as well as reconstituting lysate.
Simply add your sample, mix and incubate!
- Highly stable Control Standard Endotoxin (CSE) provided with the kit
- Cost effective premium quality reagents compared to those of competitors

INTENDED USE :

- Qualitative detection of endotoxin
- Pyrogen control or endotoxin detection during production processes

PRINCIPLE :

- Presence of the endotoxin activates factor C in the horseshoe crab amoebocyte lysate
- Activated factor C further activates factor B
- Pre-clotting enzyme (inactive) is cleaved by activated factor B into clotting enzyme (active)
- Active enzyme reduces Coagulogen into Coagulin thereby leading to clot formation

BACTERIAL ENDOTOXIN TEST (BET) Series Product Range

Product No.	Product Name	Description	Packing
BACTERIAL ENDOTOXIN TEST (BET) KITS			
THG10-0500	Endotoxin Assay Kit (gel-clot)	Amoebocyte lysate (0.5 EU/ mL) 0.2 mL x 44 pcs; BET water 2 mL x 5 pcs; Control Standard Endotoxin x 1 pcs; 44 Tests/ kit	Ampoule
THG10-0250	Endotoxin Assay Kit (gel-clot)	Amoebocyte lysate (0.25 EU/ mL) 0.2 mL x 44 pcs; BET water 2 mL x 5 pcs; Control Standard Endotoxin x 1 pcs; 44 Tests/ kit	Ampoule
THG10-0125	Endotoxin Assay Kit (gel-clot)	Amoebocyte lysate (0.125 EU/ mL) 0.2 mL x 44 pcs; BET water 2 mL x 5 pcs; Control Standard Endotoxin x 1 pcs; 44 Tests/ kit	Ampoule
THG10-0060	Endotoxin Assay Kit (gel-clot)	Amoebocyte lysate (0.06 EU/ mL) 0.2 mL x 44 pcs; BET water 2 mL x 5 pcs; Control Standard Endotoxin x 1 pcs; 44 Tests/ kit	Ampoule
THG10-0030	Endotoxin Assay Kit (gel-clot)	Amoebocyte lysate (0.03 EU/ mL) 0.2 mL x 44 pcs; BET water 2 mL x 5 pcs; Control Standard Endotoxin x 1 pcs; 44 Tests/ kit	Ampoule
HORSESHOE CRAB AMOEOCYTE LYSATE			
HG020500	Amoebocyte lysate (gel-clot)	0.5 EU/ mL, 0.2 mL x 50 pcs/ box	Ampoule
HG020250	Amoebocyte lysate (gel-clot)	0.25 EU/ mL, 0.2 mL x 50 pcs/ box	Ampoule
HG020125	Amoebocyte lysate (gel-clot)	0.125 EU/ mL, 0.2 mL x 50 pcs/ box	Ampoule
HG020060	Amoebocyte lysate (gel-clot)	0.06 EU/ mL, 0.2 mL x 50 pcs/ box	Ampoule
HG020030	Amoebocyte lysate (gel-clot)	0.03 EU/ mL, 0.2 mL x 50 pcs/ box	Ampoule
ACCESSORIES			
ET-20	Control Standard Endotoxin	10-50 EU, 50 pcs/ box	Ampoule
ET-120	Control Standard Endotoxin	51-200 EU, 50 pcs/ box	Ampoule
DCV-2000	Endotoxin indicator	2000 EU - 5000 EU, 18 pcs/ box	Ampoule
DCV50-2000	Endotoxin indicator	2000 EU - 5000 EU, 50 pcs/ box	Ampoule
DCV-4000	Endotoxin indicator	4000 EU, 18 pcs/ box	Ampoule
DCV50-4000	Endotoxin indicator	4000 EU, 50 pcs/ box	Ampoule

BACTERIAL ENDOTOXIN TEST (BET) Series Product Range

Product No.	Product Name	Description	Packing
ACCESSORIES			
WP-02	BET water	2 mL x 50 pcs/ box, endotoxin free (< 0.001 EU/ mL)	Ampoule
WP-05	BET water	5 mL x 10 pcs/ box, endotoxin free (< 0.001 EU/ mL)	Ampoule
WP-50	BET water	50 mL x 20 pcs/ box, endotoxin free (< 0.001 EU/ mL)	Vulcanized glass bottle
WP-100	BET water	100 mL x 20 pcs/ box, endotoxin free (< 0.001 EU/ mL)	Vulcanized glass bottle
PH-S-2	pH buffer (acidity)	2 mL x 50 pcs/ box, endotoxin free	Ampoule
PH-J-2	pH buffer (alkaline)	2 mL x 50 pcs/ box, endotoxin free	Ampoule
WP-SG-2 *	Pyrogen-free Ampoule	2 mL x 50 pcs/ box, endotoxin free, glucan free	-
WP-SG-5 *	Pyrogen-free Ampoule	5 mL x 10 pcs/ box, endotoxin free, glucan free	-
WP-XZ-025 *	Pyrogen-free pipette tips	250 uL x 200 pcs/ box, plastic, endotoxin free, glucan free	-
WP-XZ-100 *	Pyrogen-free pipette tips	1000 uL x 100 pcs/ box, plastic, endotoxin free, glucan free	-
GT-1285	Borosilicate glass tubes	10 x 85 mm, 40 pcs/ box	-
AE-A-2 *	AntiEnhance	2 mL x 10 pcs/ box, endotoxin free	Ampoule
For the best performance it is recommended to store above products at 2-8°C. (except codes marked with an '*')			

HiMedia Laboratories Pvt. Ltd.
www.himedialabs.com

CORPORATE OFFICE -

Plot No. C40, Road No. 21Y, MIDC, Wagle Industrial Area,
Thane (West) - 400604, Maharashtra, India.

Tel : +91-22-6147 1919 / 6116 9797 / 6903 4800 | Fax : +91-22-6147 1920

Email : info@himedialabs.com, support@himedialabs.com



For iOS



For Android



For Website