



Technical Data

Linco T Supplement

FD026

An antibiotic supplement recommended for the selective isolation and cultivation of *Neisseria* species.

Composition

Per vial sufficient for 500 ml medium

*Ingredients	Concentration
Lincomycin	0.500mg
Colistin	3mg
Amphotericin B	0.500mg
Trimethoprim lactate	3.250mg

Directions:

Rehydrate the contents of 1 vial aseptically with 10 ml of sterile distilled water. Mix gently, avoiding frothing to dissolve and aseptically add it to 240 ml of sterile, molten, cooled (45-50°C), GC Agar Base [M434](#) / GC HiVeg™ Agar Base [MV434](#) along with 250 ml of sterile 2% w/v FO Growth Supplement [FD022](#) . Also add sterile rehydrated contents of 1 vial of Yeast Autolysate Supplement [FD027](#). Mix gently and pour into sterile petri plates.

Type of specimen

Clinical samples - Stool, urine, respiratory exudates, etc.

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (1,2). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning & Precautions

In Vitro diagnostic use only. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Storage and Shelf Life

Store at 2-8°C. Use before expiry date on the label.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (1,2).

Reference

- 1.Isenberg (Ed.),2004, Clinical Microbiology Procedures Handbook, Vol.3, American Society for Microbiology, Washington. D.C.
- 2.Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology,11th Edition. Vol. 1.

* Not For Medicinal Use


Revision : 03/2022




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



In vitro diagnostic
medical device

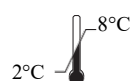

Storage temperature

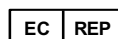

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



Do not use if
package is damaged

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal diagnostic or therapeutic use but for laboratory, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.