

LB-Growth Top Agar

G009

LB-Growth Top Agar is recommended for the cultivation and maintenance of recombinant strains of *Escherichia coli* for genetic and molecular biology studies.

Composition** :

| Ingredients | Grams/Litre |
|-----------------|-------------|
| Tryptone | 10.00 |
| Yeast extract | 5.00 |
| Sodium chloride | 10.00 |
| Agar | 7.00 |

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 32 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

Principle and Interpretation :

LB-Growth Top Agar is recommended for the cultivation and maintenance of recombinant strains of *Escherichia coli* strains for the preparation of phage and plasmid DNA. LB Media, originally developed by Bertani, is a very common and nutritionally rich growth media for *E. coli* (1). This media is extensively used for the maintenance and propagation of plasmid DNA and for the growth of the recombinant strains which contain protein expression vector (2). The ingredients provide all the growth factors required for the *E. coli* strains. Tryptone provides all the required peptides and peptones. Yeast extract supplies all the essential vitamins, nucleotides, amino acids, carbohydrate and trace elements. Sodium chloride provides the sodium ions which help in the membrane transport and maintenance of osmotic equilibrium of the medium. Top agar is used to distribute bacteriophage and bacterial cells uniformly on the thin layer over the surface of a plate. Top agar contains less amount of agar than usual plates and so stays in a molten state for several days when it is kept at 45° to 50°C.

Quality Control :

Appearance of Powder :

Cream to yellow coloured, homogeneous, free flowing powder.

Gelling :

Firm, comparable with 0.7% Agar gel.

Please refer disclaimer Overleaf

Colour and Clarity :

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates.

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Organisms (ATCC)

Escherichia coli ATCC 23724

Escherichia coli ATCC 25922

Escherichia coli MTCC1652

Growth

good-luxuriant

good-luxuriant

good-luxuriant


Storage and Shelf-life:

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Bertani, G. (1951). Studies on lysogenesis. I. The mode of phage liberation by lysogenic *Escherichia coli*. J. Bacteriol. 62:293-300.

2. Sambrook, J. E. F. Fritsch, and T. Maniatis (1989). Molecular cloning: a laboratory manual, 2nd edition ed., Cold Spring Harbour laboratory, Cold Spring Harbour, N.Y.

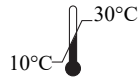

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



05/2025

PIG009_2/0522

G009-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com