

Tryptone Growth Agar II

G011

Tryptone Growth Agar II is used as a general purpose media for the growth of *Escherichia coli*.

Composition** :

Ingredients	Grams/Litre
Tryptone	10.00
Sodium chloride	8.00
Agar	15.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 33 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

Principle and Interpretation :

Tryptone Growth Agar II is used as a general purpose media for the growth of *Escherichia coli*. Tryptone provides all the required peptides and peptones. Sodium chloride provides the sodium ions which help in the membrane transport and maintenance of osmotic equilibrium of the medium.

Quality Control :

Appearance of Powder :

Cream to yellow coloured, homogeneous, free flowing powder.

Gelling :

Firm, comparable with 1.5% Agar gel.

Colour and Clarity :

Light yellow coloured, clear to slight opalescent gel forms in petriplates.

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Please refer disclaimer Overleaf

Registered Office

HiMedia Laboratories Pvt Ltd.

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,
Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797
Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920
Web : www.himedialabs.com
Email : info@himedialabs.com
mb@himedialabs.com




Organisms (ATCC)	Inoculum (CFU)	Growth
<i>Escherichia coli</i> ATCC 23724	50 - 100	good-luxuriant
<i>Escherichia coli</i> ATCC 25922	50 - 100	good-luxuriant
<i>Escherichia coli</i> MTCC 1652	50 - 100	good-luxuriant

Storage and Shelf-life :

Store below 30°C in tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Bertani, G. (1951). Studies on lysogenesis. I. The mode of phage liberation by lysogenic *Escherichia coli*. J. Bacteriol. 62:293-300.
2. Sambrook, J. E. F. Fritsch, and T. Maniatis (1989). Molecular cloning: a laboratory manual, 2nd edition ed., Cold Spring Harbour laboratory, Cold Spring Harbour, N.Y.

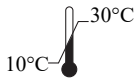

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



05/2025

PIG011_2/0522

G011-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com