

M9 Minimal Medium Salts (5X)

G013

M9 Minimal Medium Salts (5X) is recommended for growing *Escherichia coli* for molecular biology purpose.

Composition** :

Ingredients	Grams/Litre
Disodium hydrogen phosphate	33.90
Potassium dihydrogen phosphate	15.00
Sodium chloride	2.50
Ammonium chloride	5.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 56.4 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Dispense as desired and sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle and Interpretation :

M9 Minimal Salts, 5x is used in preparing M9 Minimal medium used for the cultivation and maintenance of *E. coli* in molecular biology. M9 is a chemically defined minimal growth medium for *E. coli* cultures (1, 2). It can be supplemented with specific amino acids or other required nutrients for the selection of some specific auxotroph. Sometimes the incorporation of certain additives (e.g. thiamine or casamino acids) enhances the bacterial growth. In 1959, Jacob and Monod published an article where they discussed about minimal media of *E. coli* (3). The components of this media supply the required nutrients for proper growth. Ammonium chloride acts as the nitrogen source and Sodium chloride maintains the final osmolarity of the medium. M9 Minimal Salts Base, 5x, is a 5x concentrate and should be diluted to 1x concentration before usage. Glucose may be added as a carbon source. Some other nutritional elements (e.g. calcium and magnesium) can be added as supplement.

Quality Control :

Appearance of Powder :

White to cream coloured, homogeneous, free flowing powder.

Colour and Clarity :

Colourless, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Please refer disclaimer Overleaf

Organisms (ATCC)

Escherichia coli ATCC 23724

Escherichia coli ATCC 25922

Escherichia coli MTCC1652

Growth

good-luxuriant

good-luxuriant


good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Davis, L.G., M.D. Dibner and J.F. Battey, Basic methods in molecular biology, Elsevier, new York, (1986).
2. Sambrook, J., E. F. Fritsch, and T. Maniatis, 1989, Molecular cloning: a laboratory manual, 2nd edition ed., Cold Spring Harbour laboratory, Cold Spring Harbour, N.Y.
3. Pardee, A. B., F. Jacob, and J. Monod. 1959. The genetic control and cytoplasmic expression of "inducibility" in the synthesis of β -galactosidase in *E. coli*. *J. Mol. Biol.* 1:165-178.

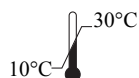

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



02/2025

PIG013_2/0222

G013-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com