

SOB Growth Medium

G014

SOB Growth Medium is used for preparing competent cells prior to transformation.

Composition** :

Ingredients	Grams/Litre
Tryptone	20.00
Yeast extract	5.00
Sodium chloride	0.50
MgSO ₄ . 7H ₂ O	5.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 30.5 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Dispense as desired and sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle and Interpretation :

SOB Growth Medium is used for preparing competent cells prior to transformation. Super Optimal Broth (SOB medium) Growth Medium was first formulated by Douglas Hanahan in 1983 as a nutritionally enriched growth media for bacteria, especially *Escherichia coli* (1). It is used in the preparation and transformation of chemically competent cells. *E. coli* cells are first made competent during transformation where perforations are made in the bacterial cells so that the foreign DNA can penetrate the cells (2). To endure this process, competent cells require a rich isotonic medium. Tryptone provides nitrogen, amino acids and other growth factors which permit the cells to go through the stress of transformation. Vitamins and trace elements are contained in Yeast Extract. Sodium chloride provides essential ions for transport and osmotic balance. Magnesium sulfate provides magnesium ions which are required in a variety of enzymatic reactions, including DNA replication (2).

Quality Control :

Appearance of Powder :

Light yellow coloured, homogeneous, free flowing powder.

Colour and Clarity :

Light amber coloured, clear solution without any precipitate.

Please refer disclaimer Overleaf

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Organisms (ATCC)

Escherichia coli ATCC 23724

Escherichia coli ATCC 25922

Escherichia coli MTCC1652

Growth

good-luxuriant

good-luxuriant


good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Hanahan, D., 1983. J. Mol. Biol.166:557.
2. Sambrook J., E. F. Fritsch, and T. Maniatis. 1989. Molecular cloning: a laboratory manual, 2nd ed. Cold Spring Harbor Laboratory, Cold Spring Harbor, N.Y.

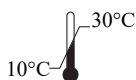

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



05/2025

PIG014_2/0522

G014-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com