

NZYDT Growth Medium

G026

NZYDT Growth Medium is used for cultivation of recombinant *Escherichia coli* strains.

Composition** :

Ingredients	Grams/Litre
Casein enzymic hydrolysate	10.00
Yeast extract	5.00
MgSO ₄ . 7H ₂ O	1.00
Sodium chloride	5.00
Thymidine	0.04

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 20.53 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Dispense as desired and sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle and Interpretation :

NZYDT Growth Medium is used for cultivation of recombinant *Escherichia coli* strains. This medium was developed by Blattner and colleagues as a rich medium for the propagation of bacteriophages which infects its cell host, *E. coli* for replication (1). Cells grow very fast in this medium as this medium provides all the amino acids, vitamins and other metabolites required for cell growth (2). Casein enzymic hydrolysate provides nitrogen, amino acids, and carbon sources for the cells. Yeast extract functions as the source of vitamins and trace elements. Sodium chloride provides sodium ions for transport and osmotic balance and Magnesium sulfate is a source of magnesium ions required in a variety of enzymatic reactions, including DNA replication (3). NZYDT Growth Medium allows the cells to grow more rapidly as they do not have to synthesize nucleotide precursors and other factors required for growth.

Quality Control :

Appearance of Powder :

Light yellow coloured, homogeneous, free flowing powder.

Please refer disclaimer Overleaf

1



Registered Office

HiMedia Laboratories Pvt Ltd.

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,
Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797
Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920

Web : www.himedialabs.com

Email : info@himedialabs.com

mb@himedialabs.com

The information contained herein is believed to be accurate and complete. However no warranty or guarantee whatsoever is made or is to be implied with respect to such information or with respect to any product, method or apparatus referred to herein

Colour and Clarity :

Light yellow coloured, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Organisms (ATCC)

Escherichia coli ATCC 23724

Escherichia coli ATCC 25922

Escherichia coli MTCC1652

Growth

good-luxuriant

good-luxuriant


good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

- (1) Blattner, F. R., B. G. Williams, A. E. Blechl, K. Denniston-Thompson, H. E. Faber, L. A. Furlong, D. J. Grunwald, D. O. Kiefer, D. D. Moore, J. W. Schumm, E. L. Sheldon, and O. Smithies. 1977. Charon phages: Safer derivatives of bacteriophage for DNA cloning. *Science* 196:161.
- (2) Ausubel, F. M., R. Brent, R. E. Kingston, D. D. Moore, J. G. Seidman, J. A. Smith, and K. Struhl. 1994. *Current protocols in molecular biology*, vol. 1. Current Protocols, New York, NY.
- (3) Sambrook J., E. F. Fritsch, and T. Maniatis. 1989. *Molecular cloning: a laboratory manual*, 2nd ed. Cold Spring Harbor Laboratory, Cold Spring Harbor, NY.

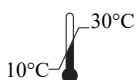

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
 Reg. Off: Plot No. C-40, Road No. 21Y,
 MIDC, Wagle Industrial Area, Thane,
 (West) 400604, Maharashtra, INDIA.
 Web: www.himedialabs.com



04/2027

PIG026_2/0424

G026-03

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg. office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
 Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com