

YPD (YEPD) Growth Agar

G038

YPD (YEPD) Growth Agar is used for the growth of *Saccharomyces cerevisiae*.

Composition** :

Ingredients	Grams/Litre
Peptone	20.00
Yeast extract	10.00
Dextrose	20.00
Agar	15.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 65 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

Principle and Interpretation :

YPD (YEPD) Growth Agar is used for the growth of *Saccharomyces cerevisiae*. Yeasts are unicellular eukaryotes and extensively studied model organism in molecular genetics. They are chemoorganotrophs as they utilize organic compounds as a source of energy. YPD (YEPD) Growth Agar is used for the maintenance and propagation of yeasts including *S. cerevisiae* in various molecular microbiology procedures (1, 2). YPD functions as a complete medium for yeast growth and it contains yeast extract, peptone and glucose or dextrose. Yeast extract supplies B-complex vitamins and it contains all the amino acids necessary for growth. Peptone acts as the source of nitrogen, vitamins and minerals. Dextrose serves as the carbon source. This medium supports the vigorous growth of wild type as well as mutant strains of all kinds of budding yeast.

Quality Control :

Appearance of Powder :

Light yellow coloured, homogeneous, free flowing powder.

Gelling :

Firm, comparable with 1.5% Agar gel

Colour and Clarity :

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates.

Please refer disclaimer Overleaf

1



Registered Office

HiMedia Laboratories Pvt Ltd.

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,
Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797
Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920

Web : www.himedialabs.com
Email : info@himedialabs.com
mb@himedialabs.com

The information contained herein is believed to be accurate and complete. However no warranty or guarantee whatsoever is made or is to be implied with respect to such information or with respect to any product, method or apparatus referred to herein

Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Organisms (ATCC)

Saccharomyces cerevisiae ATCC9763

Growth


good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Adams, A., D. E. Gottschling, C. A. Kaiser, and T. Stearns. 1997. Methods in yeast genetics: A Cold Spring Harbor Laboratory Course Manual. Cold Spring Harbor Laboratory Press, Cold Spring Harbor, New York.
2. Burke, D., Dawson, D., and T. Stearns. 2000. Method in yeast genetics. Cold Spring Harbor Laboratory Press, Cold Spring Harbor, New York.

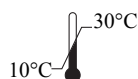

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



12/2024

PIG038_1/1221

G038-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com