

## Sporulation Growth Agar

G042

Sporulation Growth Agar is used for the growth and sporulation of *Saccharomyces cerevisiae*.

### Composition\*\* :

Ingredients	Grams/Litre
Yeast extract	1.00
Potassium acetate	20.00
Agar	15.00

\*\* Formula adjusted, standardized to suit performance parameters

### Directions :

Suspend 36 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

### Principle and Interpretation :

Sporulation Growth Agar is used for the growth and sporulation of *Saccharomyces cerevisiae*. *Saccharomyces cerevisiae* is a unicellular eukaryote that has become an important tool in microbial genetic techniques. It undergoes meiosis and sporulation which takes place in a single cell. This medium stimulates the sporulation of diploid yeast cells which occurs only when cells are deprived of any carbon source. Potassium acetate enhances the sporulation of diploid strains. This medium contains agar as the solidifying agent.

### Quality Control :

#### Appearance of Powder :

Cream to yellow coloured, homogeneous, free flowing powder.

#### Gelling :

Firm, comparable with 1.5% Agar gel.

#### Colour and Clarity :

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates.

#### Cultural Response :

Cultural characteristics observed after an incubation at 35-37°C for 18 - 48 hours.

Please refer disclaimer Overleaf

#### Registered Office

#### HiMedia Laboratories Pvt Ltd.

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,  
Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797  
Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920  
Web : [www.himedialabs.com](http://www.himedialabs.com)  
Email : [info@himedialabs.com](mailto:info@himedialabs.com)  
[mb@himedialabs.com](mailto:mb@himedialabs.com)



**Organisms (ATCC)**  
*Saccharomyces cerevisiae* ATCC9763


**Growth**  
good-luxuriant

**Storage and Shelf-life :**

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

**References:**

1. Adams, A., D. E. Gottschling, C. A. Kaiser, and T. Stearns. 1997. Methods in yeast genetics: A Cold Spring Harbor Laboratory Course Manual. Cold Spring Harbor Laboratory Press, Cold Spring Harbor, New York.
2. Burke, D., Dawson, D., and T. Stearns. 2000. Method in yeast genetics. Cold Spring Harbor Laboratory Press, Cold Spring Harbor, New York.

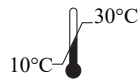

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: [www.himedialabs.com](http://www.himedialabs.com)



05/2025

PIG042\_1/0522

G042-02

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: [techhelp@himedialabs.com](mailto:techhelp@himedialabs.com) Website: [www.himedialabs.com](http://www.himedialabs.com)