

EMM Growth Medium with Sorbitol

G056

EMM Growth Medium with Sorbitol is a minimal defined media for the growth of *Schizosaccharomyces pombe*.

Composition :**

Ingredients	Grams/Litre
Phthalic acid K+	3.00
Disodium hydrogen phosphate	2.20
Ammonium chloride	5.00
Sorbitol	216.00
Magnesium chloride, 6H ₂ O	1.05
Calcium chloride, 2H ₂ O	0.0147
Potassium chloride	1.00
Sodium sulphate	0.04
Pantothenic acid	0.001
Nicotinic acid	0.01
Myoinositol	0.01
Biotin	0.001
Boric acid	0.0005
Manganese sulphate	0.0004
Zinc sulphate, 7H ₂ O	0.0004
Ferric chloride, 6H ₂ O	0.0002
Molybdic acid	40 mcg
Potassium iodide	0.0001
Copper sulphate, 5H ₂ O	40 mcg
Citric acid	0.001

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 228.35 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

Principle and Interpretation :

EMM Growth Medium with Sorbitol is a minimal defined media for the growth of *Schizosaccharomyces pombe*. Yeasts are unicellular eukaryotes and extensively studied model organism in molecular genetics. The fission yeast *Schizosaccharomyces pombe* is a model eukaryote which is very useful in studies of cell cycle and chromosome dynamics. These cells maintain their shape by growing through the cell tips and divide by medial fission to produce two daughter cells of equal sizes that makes them a powerful tool in cell cycle research. It was first developed as an experimental model in the 1950's for studying genetics (1, 2) and for studying the cell cycle (3, 4). EMM (Edinburgh Minimal Media) Growth Medium with Sorbitol is used for the maintenance and propagation of *S. pombe* in various molecular microbiology procedures. It functions as a minimal defined medium for fission yeast growth and it contains sorbitol, minerals and trace elements.

Quality Control :**Appearance of Powder :**

Cream to yellow coloured, homogeneous, free flowing powder.

Colour and Clarity :

Light yellow coloured, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

Organisms

Schizosaccharomyces pombe

Growth


good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Leupold U. (1950) CR Trav Lab Carlsberg Ser Physiol 24:381-480.
2. Leupold U. (1993) The origins of *Schizosaccharomyces pombe* genetics. In: Hall MN, Linder P. eds. The early Days of Yeast Genetics. New York. Cold Spring Harbor Laboratory Press. 125-128.
3. Mitchinson JM. (1975) Exp Cell Res 13:244-262.
4. Mitchinson JM. (1990) Bioessays 4:189-191.

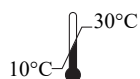

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



08/2025

PIG056_2/0822

G056-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com