

## EMM Growth Medium with Low Dextrose

G057

EMM Growth Medium with Low Dextrose is a minimal defined media for the growth of *Schizosaccharomyces pombe*.

### Composition\*\* :

Ingredients	Grams/Litre
Phthalic acid K+	3.00
Disodium hydrogen phosphate	2.20
Ammonium chloride	5.00
Dextrose	5.00
Magnesium chloride, 6H <sub>2</sub> O	1.05
Calcium chloride, 2H <sub>2</sub> O	0.0147
Potassium chloride	1.00
Sodium sulphate	0.04
Pantothenic acid	0.001
Nicotinic acid	0.01
Myoinositol	0.01
Biotin	0.001
Boric acid	0.0005
Manganese sulphate	0.0004
Zinc sulphate, 7H <sub>2</sub> O	0.0004
Ferric chloride, 6H <sub>2</sub> O	0.0002
Molybdcic acid	40 mcg
Potassium iodide	0.0001
Copper sulphate, 5H <sub>2</sub> O	40 mcg
Citric acid	0.001

\*\* Formula adjusted, standardized to suit performance parameters

### Directions :

Suspend 17.2 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

### Principle and Interpretation :

EMM Growth Medium with Low Dextrose is a minimal defined media for the growth of *Schizosaccharomyces pombe*.

---

Please refer disclaimer Overleaf

1



Registered Office

**HiMedia Laboratories Pvt Ltd.**

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,  
Thane, (West) 400604, Maharashtra, INDIA.

Customer Care No.: 00-91-22-6116 9797

Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920

Web : [www.himedialabs.com](http://www.himedialabs.com)

Email : [info@himedialabs.com](mailto:info@himedialabs.com)

[mb@himedialabs.com](mailto:mb@himedialabs.com)

The information contained herein is believed to be accurate and complete. However no warranty or guarantee whatsoever is made or is to be implied with respect to such information or with respect to any product, method or apparatus referred to herein

**Quality Control :**

**Appearance of Powder :**

Light yellow coloured, homogeneous, free flowing powder.

**Colour and Clarity :**

Light yellow coloured, clear solution without any precipitate.

**Cultural Response :**

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

**Organisms (ATCC)**

*Schizosaccharomyces pombe*

**Growth**


good-luxuriant

**Storage and Shelf-life:**

Upon receipt, store at 2 - 8°C. Use before expiry date on the label.

**References:**

1. Leupold U. (1950) CR Trav Lab Carlsberg Ser Physiol 24:381-480.
2. Leupold U. (1993) The origins of Schizosaccharomyces pombe genetics. In: Hall MN, Linder P. eds. The early Days of Yeast Genetics. New York. Cold Spring Harbor Laboratory Press. 125-128.
3. Mitchinson JM. (1975) Exp Cell Res 13:244-262.
4. Mitchinson JM. (1990) Bioessays 4:189-191.

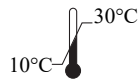

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: www.himedialabs.com



04/2025

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.