

**SD Growth Medium**

**G062**

SD Growth Medium is a synthetic defined media for the growth of *Saccharomyces cerevisiae*.

**Composition\*\* :**

<b>Ingredients</b>	<b>Grams/Litre</b>
Potassium dihydrogen phosphate	1.00
Magnesium sulphate	0.50
Sodium chloride	0.10
Calcium chloride	0.10
Biotin	0.002 mg
Calcium pantothenate	0.4 mg
Folic acid	0.002 mg
Inositol	2.00 mg
Niacin	0.4 mg
PABA	0.2 mg
Pyridoxin, HCl	0.4 mg
Riboflavin	0.2 mg
Thiamine HCl	0.4 mg
Boric acid	0.5 mg
Copper sulphate	0.04 mg
Potassium iodide	0.1 mg
Ferric chloride	0.2 mg
Manganese sulphate	0.4 mg
Sodium molybdate	0.2 mg
Zinc sulphate	0.4 mg
Ammonium sulphate	5.00
Dextrose	20.00
Adenine	0.010
L-Arginine HCl	0.050
L-Aspartic acid	0.080
L-Histidine HCl	0.020
L-Isoleucine	0.050
L-Leucine	0.100
L-Lysine HCl	0.050
L-Methionine	0.020
L-Phenylalanine	0.050
L-Threonine	0.100
L-Tryptophan	0.050
L-Tyrosine	0.050

## SD Growth Medium

G062

Uracil	0.020
L-Valine	0.140

\*\* Formula adjusted, standardized to suit performance parameters

### Directions :

Suspend 27.49 grams in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

### Principle and Interpretation :

SD Growth Medium is a synthetic defined agar media for the growth of *Saccharomyces cerevisiae*. Synthetically Defined media known as Yeast Nitrogen Base Media for the growth of yeast cells were first cited by Wickerham (1, 2). SD Growth Medium includes a yeast nitrogen base along with ammonium sulfate and dextrose as the carbon source, which is further supplemented with various nutrients. This makes it a complete synthetic growth medium for yeast cells.

### Quality Control :

#### Appearance of Powder :

Off white to yellow coloured, homogeneous, free flowing powder.

#### Gelling :

Firm, comparable with 1.5% Agar gel.

#### Colour and Clarity :

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates.

#### Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

#### Organisms (ATCC)

*Saccharomyces cerevisiae* ATCC 9763

#### Growth

good-luxuriant

---

Please refer disclaimer Overleaf

## SD Growth Medium


G062

### Storage and Shelf-life :

Upon receipt, store at 2 - 8°C. Use before expiry date on the label.

### References :

1. Wickerham L. J., 1951, U.S. Dept. Agric. Tech. Bull. No. 1029
2. Wickerham L. J., 1946, J. Bacteriol., 52:293
3. Kaiser, C., et al. Methods in Yeast Genetics Cold Spring Harbor, (1994)

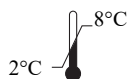

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: [www.himedialabs.com](http://www.himedialabs.com)



04/2025

PIG062\_2/0422

G062-02

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: [techhelp@himedialabs.com](mailto:techhelp@himedialabs.com) Website: [www.himedialabs.com](http://www.himedialabs.com)