

**SD Growth Medium w/o HIS-LEU-TRP**

**G070**

SD Growth Medium w/o HIS-LEU-TRP is a synthetic defined media for the selective growth of *Saccharomyces cerevisiae*.

**Composition\*\* :**

<b>Ingredients</b>	<b>Grams/Litre</b>
Potassium dihydrogen phosphate	1.00
Magnesium sulphate	0.50
Sodium chloride	0.10
Calcium chloride	0.10
Biotin	0.002 mg
Calcium pantothenate	0.4 mg
Folic acid	0.002 mg
Inositol	2.00 mg
Niacin	0.4 mg
PABA	0.2 mg
Pyridoxin, HCl	0.4 mg
Riboflavin	0.2 mg
Thiamine HCl	0.4 mg
Boric acid	0.5 mg
Copper sulphate	0.04 mg
Potassium iodide	0.1 mg
Ferric chloride	0.2 mg
Manganese sulphate	0.4 mg
Sodium molybdate	0.2 mg
Zinc sulphate	0.4 mg
Ammonium sulphate	5.00
Dextrose	20.00
Adenine	0.010
L-Arginine HCl	0.050
L-Aspartic acid	0.080
L-Isoleucine	0.050
L-Lysine HCl	0.050
L-Methionine	0.020
L-Phenylalanine	0.050
L-Threonine	0.100
L-Tyrosine	0.050
Uracil	0.020
L-Valine	0.140

\*\* Formula adjusted, standardized to suit performance parameters

**Directions :**

Suspend 27.32 grams in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

**Principle and Interpretation :**

SD Growth Medium w/o HIS-LEU-TRP is a synthetic defined media for the selective growth of *Saccharomyces cerevisiae*. Synthetically Defined media known as Yeast Nitrogen Base Media for the growth of yeast cells were first cited by Wickerham (1, 2). Synthetic defined Growth Medium w/o HIS-LEU-TRP includes a yeast nitrogen base along with ammonium sulfate and dextrose as the carbon source, which is further supplemented with various amino acids except histidine, leucine and tryptophan. This makes it a triple dropout growth medium for yeast cells. A histidine, leucine and tryptophan auxotrophic yeast mutant strain cannot grow on this media but a wild-type or a histidine, leucine and tryptophan prototrophic yeast strain can grow. The histidine, leucine and tryptophan auxotroph has mutation in the genes (e.g. *HIS3*, *LEU2* and *TRP1*) of the histidine, leucine as well as tryptophan synthesis pathway and this mutant strain will grow in this medium if histidine, leucine and tryptophan are supplied from outside e.g. from plasmids which contain *HIS3*, *LEU2* and *TRP1* gene (3). For this purpose, a *his3leu2trp1* mutant strain of *S. cerevisiae* is transformed with *HIS3*, *LEU2* and *TRP1* containing plasmid and the transformants can be selected by growing the cells on SD Growth Medium w/o HIS-LEU-TRP.

**Quality Control :****Appearance of Powder :**

White to cream coloured, homogeneous, free flowing powder.

## SD Growth Medium w/o HIS-LEU-TRP

G070

### Colour and Clarity :

Colourless to light yellow coloured, clear solution without any precipitate.

### Cultural Response :

Cultural characteristics observed after incubation at 25-30°C for 18 - 48 hours.

### Organisms (ATCC)

*Saccharomyces cerevisiae* ATCC9763

### Growth


good-luxuriant

### Storage and Shelf-life :

Upon receipt, store at 2 - 8°C. Use before expiry date on the label.

### References :

1. Wickerham L. J., 1951, U.S. Dept. Agric. Tech. Bull. No. 1029
2. Wickerham L. J., 1946, J. Bacteriol., 52:293
3. Kaiser, C., et al. Methods in Yeast Genetics Cold Spring Harbor, (1994)

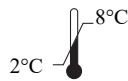

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: [www.himedialabs.com](http://www.himedialabs.com)



04/2025

PIG070\_2/0422

G070-02

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: [techhelp@himedialabs.com](mailto:techhelp@himedialabs.com) Website: [www.himedialabs.com](http://www.himedialabs.com)