

DOB Growth Medium

G082

DOB Growth Medium is used for the growth of *Saccharomyces cerevisiae*.

Composition** :

Ingredients	Grams/Litre
Potassium dihydrogen phosphate	1.00
Magnesium sulphate	0.50
Sodium chloride	0.10
Calcium chloride	0.10
Biotin	0.002 mg
Calcium pantothenate	0.4 mg
Folic acid	0.002 mg
Inositol	2.00 mg
Niacin	0.4 mg
PABA	0.2 mg
Pyridoxin, HCl	0.4 mg
Riboflavin	0.2 mg
Thiamine HCl	0.4 mg
Boric acid	0.5 mg
Copper sulphate	0.04 mg
Potassium iodide	0.1 mg
Ferric chloride	0.2 mg
Manganese sulphate	0.4 mg
Sodium molybdate	0.2 mg
Zinc sulphate	0.4 mg
Ammonium sulphate	5.00
Dextrose	20.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 26.7 grams in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

Principle and Interpretation :

DOB Growth Medium is used for the growth of *Saccharomyces cerevisiae*. This yeast strain is called budding yeast and is extensively studied microorganism in molecular and cell biology.

Please refer disclaimer Overleaf

DOB Growth Medium is developed according to the formulae of Wickerham (1, 2) which supplies the required nutritional elements for budding yeast cells. This media contains all the essential vitamins and inorganic salts which are required for the propagation of yeast cells but does not contain the amino acids namely, histidine, methionine, leucine, lysine and tryptophan. Dextrose serves as the carbon source for the cells.

Quality Control :

Appearance of Powder :

Light yellow coloured, homogeneous, free flowing powder.

Colour and Clarity :

Light yellow coloured, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

Organisms (ATCC)

Saccharomyces cerevisiae

Growth


good-luxuriant

References :

1. Wickerham. 1939. J. Tropical Med. Hyg. 42:176.
2. Wickerham. 1946. J. Bacteriol. 52:293.

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

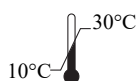

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



08/2025

PIG082_2/0822

G082-02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com