

Yeast Nitrogen Base (YNB) w/o Ammonium Sulphate

G090

Yeast Nitrogen Base (YNB) w/o Ammonium Sulphate is used for the growth of all strains of *Saccharomyces cerevisiae*.

Composition** :

Ingredients	Milligrams/Litre
Potassium phosphate, monobasic	1000.00
Magnesium sulphate	500.00
Sodium chloride	100.00
Calcium chloride	100.00
Biotin	0.002
Calcium pantothenate	0.40
Folic acid	0.002
Inositol	2.00
Niacin	0.40
PABA	0.20
Pyridoxin, HCl	0.40
Riboflavin	0.20
Thiamine HCl	0.40
Boric acid	0.50
Copper sulphate	0.04
Potassium iodide	0.10
Ferric chloride	0.20
Manganese sulphate	0.40
Sodium molybdate	0.20
Zinc sulphate	0.40

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 1.7 grams in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

Principle and Interpretation :

Yeast Nitrogen Base (YNB) w/o Ammonium Sulphate is used for the growth of all strains of *Saccharomyces cerevisiae*. This yeast strain is called budding yeast and is extensively studied microorganism in molecular and cell biology.

Yeast Nitrogen Base (YNB) w/o Ammonium Sulfate is developed according to the formulae of Wickerham (1, 2) which supplies the required nutritional elements for

Please refer disclaimer Overleaf

budding yeast cells. This media is used for the classification of yeast strains based on the carbon and nitrogen requirements and it contains all the essential vitamins and inorganic salts except Ammonium sulphate which is required for the propagation of yeast cells. Furthermore, this media does not contain the amino acids namely, histidine, methionine, leucine, lysine and tryptophan.

Quality Control :

Appearance of Powder :

Light yellow coloured, homogeneous, free flowing powder.

Colour and Clarity :

Light yellow coloured, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

Organisms (ATCC)

Saccharomyces cerevisiae

Growth


good-luxuriant

References :

1. Wickerham. 1939. J. Tropical Med. Hyg. 42:176.
2. Wickerham. 1946. J. Bacteriol. 52:293.

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

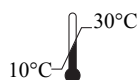

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



12/2024

PIG090_1/1221

G090_02

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com