

## Complete Supplement Mixture (CSM)

G100

Complete Supplement Mixture (CSM) is used as a dropout supplement for all strains of *Saccharomyces cerevisiae*.

### Composition\*\* :

Ingredients	Milligrams/Litre
Adenine	10.00
L-Arginine HCl	50.00
L-Aspartic acid	80.00
L-Histidine HCl	20.00
L-Isoleucine	50.00
L-Leucine	100.00
L-Lysine HCl	50.00
L-Methionine	20.00
L-Phenylalanine	50.00
L-Threonine	100.00
L-Tryptophan	50.00
L-Tyrosine	50.00
Uracil	20.00
L-Valine	140.00

\*\* Formula adjusted, standardized to suit performance parameters

### Directions :

Suspend 790 mg in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

### Principle and Interpretation :

Complete Supplement Mixture (CSM) is used as a dropout supplement for all strains of *Saccharomyces cerevisiae*. This yeast strain is called budding yeast and is extensively studied microorganism in molecular and cell biology.

Complete Supplement Mixture (CSM) is developed according to the formulae of Wickerham (1, 2) which supplies the essential nutritional elements for budding yeast cells. It is composed of all the amino acids required for the vigorous growth of *Saccharomyces cerevisiae*. It can be mixed with yeast nitrogen base (YNB), ammonium sulphate and an appropriate carbon source (glucose/galactose/raffinose) to produce a media suitable for the growth of budding yeast. It is also the basis for many “Drop-Out” formulations for the selections of auxotrophic requirements and transformants.

Please refer disclaimer Overleaf

1



Registered Office

HiMedia Laboratories Pvt Ltd.

Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area,  
Thane, (West) 400604, Maharashtra, INDIA.

Customer Care No.: 00-91-22-6116 9797

Tel : 00-91-22-6147 1919, 6903 4800

Fax : 6147 1920

Web : [www.himedialabs.com](http://www.himedialabs.com)

Email : [info@himedialabs.com](mailto:info@himedialabs.com)

[mb@himedialabs.com](mailto:mb@himedialabs.com)

## Quality Control :

### Appearance of Powder :

White to cream coloured, homogeneous, free flowing powder.

### Colour and Clarity :

Colourless, clear solution without any precipitate.

### Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

### Organisms (ATCC)

*Saccharomyces cerevisiae* ATCC 9763

### Growth


good-luxuriant

## Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

## References :

1. Wickerham. 1939. J. Tropical Med. Hyg. 42:176.
2. Wickerham. 1946. J. Bacteriol. 52:293.

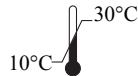

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: [www.himedialabs.com](http://www.himedialabs.com)



06/2025

PIG100\_2/0622

G100-02

## Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: [techhelp@himedialabs.com](mailto:techhelp@himedialabs.com) Website: [www.himedialabs.com](http://www.himedialabs.com)