

4 XYT Growth Medium

G151

4 XYT Growth Medium is an optimized formulation for the cultivation of recombinant *Escherichia coli*.

Composition** :

Ingredients	Grams/Litre
Tryptone	32.00
Yeast extract	20.00
Sodium chloride	5.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 57 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well and dispense as desired.

Principle and Interpretation :

4 XYT Growth Medium is an optimized formulation for the cultivation of recombinant *Escherichia coli*. This media is 4 times richer than the YT media. This media was originally formulated as a nutritionally enriched growth medium for growth of recombinant strains of *Escherichia coli* and can also be used for propagation of M13 bacteriophage (1-3). It permits larger quantity of phage production without exhausting the host. Yeast extract and tryptone provide all the required amino acids, nucleotide precursors, vitamins and other metabolites and as a result the cells grow faster in this medium. Sodium chloride provides sodium ions for transport and osmotic balance.

Quality Control :

Appearance of Powder :

Cream to yellow coloured, homogeneous, free flowing powder.

Colour and Clarity :

Light yellow to yellow coloured, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after incubation at 35-37°C for 18 - 48 hours.


Organisms (ATCC)	Growth
<i>Escherichia coli</i> ATCC 23724	good-luxuriant
<i>Escherichia coli</i> ATCC 25922	good-luxuriant
<i>Escherichia coli</i> DH5alpha MTCC1652	good-luxuriant

Storage and Shelf-life :

Store below 30°C in tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label.

References:

1. Difco manual 11th ed., Sparks, MD (1998), 22-23
2. Assubel, F.M., R. Brent, R.E. Kingston, D.D. Moore, J.G. Seidman, J.A. Smith and K. Struhl, Current protocols in molecular biology, vol. 1, Current Protocols, New York, (1994)
3. Davis, L.G., M.D. Dibner and J.F. Battey, Basic methods in molecular biology, Elsevier, new York, (1986).

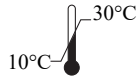

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,
Reg. Off: Plot No. C-40, Road No. 21Y,
MIDC, Wagle Industrial Area, Thane,
(West) 400604, Maharashtra, INDIA.
Web: www.himedialabs.com



04/2027

PIG151_2/0424

G151-03

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: techhelp@himedialabs.com Website: www.himedialabs.com