

Fluid Thioglycollate medium, Granulated[®]

GM009

Intended use

Recommended for sterility testing of biologicals and for cultivation of anaerobes, aerobes and microaerophiles from pharmaceutical and clinical samples.

Composition**

Ingredients	g / L
Tryptone	15.000
Yeast extract	5.000
Dextrose (Glucose)	5.500
Sodium chloride	2.500
L-Cystine	0.500
Sodium thioglycollate	0.500
Resazurin sodium	0.001
Agar	0.750
Final pH (at 25°C)	7.1±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 29.75 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 25°C and store in a cool dark place preferably below 25°C. Note : If more than the upper one-third of the medium has acquired a pink-purple colour, the medium may be restored once by heating in a water bath or in free flowing steam until the pink-purple colour disappears.

Principle And Interpretation

Brewer (1) formulated Fluid Thioglycollate Medium for rapid cultivation of aerobes as well as anaerobes including microaerophiles by adding a reducing agent and small amount of agar. The BP (2), EP (3), USP (4), and AOAC (5) have recommended the media for sterility testing of antibiotics, biologicals and foods and for determining the phenol coefficient and sporicidal effect of disinfectants. However, it is intended for the examination of clear liquid or water-soluble materials. Fluid Thioglycollate Medium is also routinely used to check the sterility of stored blood in blood banks (6). Dextrose, tryptone, yeast extract, L-cystine provide the growth factors necessary for bacterial multiplication. L-cystine and sodium thioglycollate allows *Clostridium* to grow in this medium even under aerobic conditions (7). Also the small amount of agar used in the medium favors the growth of aerobes as well as anaerobes in the medium, even if sodium thioglycollate is deleted from the medium(1). Sodium thioglycollate act as a reducing agent and neutralizes the toxic effects of mercurial preservatives and peroxides formed in the medium, thereby promoting anaerobiosis, and making the medium suitable to test materials containing heavy metals. (8,9). Any increase in the oxygen content is indicated by a colour change of redox indicator, resazurin to red (6,10,11). The small amount of agar helps in maintaining low redox potential for stabilizing the medium (9).

Type of specimen

Pharmaceutical samples for sterility testing, clinical samples- pus, wounds

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (12,13).

For pharmaceutical samples, follow appropriate techniques for sample collection, processing as per guidelines (2,3,4) After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. It is intended for the examination of clear liquid or water-soluble materials.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control**Appearance**

Cream to yellow granulated free flowing powder

Colour and Clarity of prepared medium

Light straw coloured, clear to slightly opalescent solution with upper 10% or less medium pink-purple on standing.

Reaction

Reaction of 2.97% w/v aqueous solution at 25°C. pH : 7.1±0.2

pH

6.90-7.30

Growth Promotion Test

As per USP/EP/BP/IP

Growth promoting properties

Clearly visible growth of microorganism comparable to that previously obtained with previously tested and approved lot of medium occurs at the specified temperature for not more than the shortest period of time specified inoculating ≤100 cfu at 30-35°C for or not more than 3 days for aerobes and anaerobes.

Sterility Testing + Validation

The medium is tested with suitable strains of microorganisms inoculating ≤100cfu and incubating at 20-25°C for not more than 3 days in case of bacteria and not more than 5 days in case of fungi.

Testing in accordance with EN ISO 11133:2014/Amd.1:2018(E) (10)

Cultural characteristics observed after an incubation at 36-38°C for 18-24 hours

Organism	Inoculum (CFU)	Growth	Incubation at
Growth promoting			
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 6538 (00032*)	50 -100	luxuriant	30-35°C
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*) ^	50 -100	luxuriant	30-35°C
\$ <i>Bacillus spizizenii</i> ATCC 6633 (00003*)	50 -100	luxuriant	30-35°C
^ <i>Pseudomonas paraeruginosa</i> ATCC 9027 (00026*)	50 -100	luxuriant	30-35°C
## <i>Kocuria rhizophila</i> ATCC 9341	50 -100	luxuriant	30-35°C
<i>Clostridium sporogenes</i> ATCC 19404 (00008*)	50 -100	luxuriant	30-35°C
<i>Clostridium sporogenes</i> ATCC 11437	50 -100	luxuriant	30-35°C
\$\$ <i>Phocaeicola vulgatus</i> ATCC 8482	50 -100	luxuriant	30-35°C
<i>Escherichia coli</i> ATCC 8739 (00012*)	50 -100	luxuriant	30-35°C
<i>Salmonella Typhimurium</i> ATCC 14028 (00031*)	50 -100	luxuriant	30-35°C
<i>Salmonella Abony</i> NCTC 6017	50 -100	luxuriant	30-35°C
Sterility Testing- Growth promotion + Validation			
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 6538 (00032*)	50 -100	luxuriant	20-25°C
\$ <i>Bacillus spizizenii</i> ATCC 6633 (00003*)	50 -100	luxuriant	20-25°C

[^] <i>Pseudomonas paraeruginosa</i> ATCC 9027 (00026*)	50 -100	luxuriant	20-25°C
^{##} <i>Kocuria rhizophila</i> ATCC 9341	50 -100	luxuriant	20-25°C
<i>Candida albicans</i> ATCC 10231 (00054*)	50 -100	luxuriant	20-25°C
[#] <i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	50 -100	luxuriant	20-25°C

Testing in accordance with EN ISO 11133:2014/Amd.1:2018(E) (10)

Cultural characteristics observed after an incubation at 36-38°C for 18-24 hours

<i>Clostridium perfringens</i> ATCC 13124 (00007*)	50 -100	luxuriant	36-38°C
---	---------	-----------	---------

Key : * Corresponding WDCM numbers, # Formerly known as *Aspergillus niger*,
 ## Formerly known as *Micrococcus luteus* \$ Formerly known as *Bacillus subtilis* subsp. *spizizenii*
 \$\$ Formerly known as *Bacteroides vulgatus* ^ Formerly known as *Pseudomonas aeruginosa*

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (12,13).

Reference


1. Brewer, 1940, J. Am. Med. Assoc., 115:598.
2. The British Pharmacopoeia, 2022, Medicines and Healthcare products Regulatory Agency.
3. European Pharmacopoeia, 2022, 10 th volume, European Directorate for the quality of medicines & Healthcare.
4. The United States Pharmacopoeia-National Formulary (USP-NF), 2022
5. Williams H., (Ed.), 2005, Official Methods of Analysis of the Association of Official Analytical Chemists, 19th Ed., AOAC, Washington, D.C
6. Nungester, Hood and Warren, 1943, Proc. Soc. Exp. Biol. Med., 52:287.
7. Portwood, 1944, J. Bact., 48:255.
8. Federal Register, 1992, Fed. Regist., 21:640.
9. Quastel and Stephenson, 1926, J.Biochem., 20
10. Marshall, Gunnison and Luxen, 1940, Proc. Soc. Exp. Biol. Med., 43:672.
11. MacFaddin J.F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams and Wilkins, Baltimore.
12. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
13. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC,WagleIndustrial Area,
Thane (W) -400604, MS, India



IVD *In vitro diagnostic
medical device*

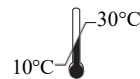

Storage temperature

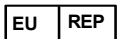

AR Experts BV
Boeingavenue 209
1119 PD Schiphol-Rijk
The Netherlands



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.