



Dichloran Glycerol Medium Base, Granulated[®]

GM1129

Intended use

Recommended for selective isolation of xerophilic moulds from food and clinical samples.

Composition**

Ingredients	g / L
Peptone	5.000
Dextrose (Glucose)	10.000
Potassium dihydrogen phosphate	1.000
Magnesium sulphate	0.500
Dichloran	0.002
Chloramphenicol	0.100
Agar	15.000
Final pH (at 25°C)	5.6±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 15.8 grams in 500 ml purified / distilled water. Heat to boiling to dissolve the medium completely. Add 110 grams of glycerol (Analytical Reagent Grade). Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

Principle And Interpretation

Dichloran Glycerol Medium was formulated by Hocking and Pitt (1) and is recommended for isolation and enumeration of xerophilic moulds from dried and semidried foods. The glycerol at 18% (w/w) lowers the water activity (aw) from 0.999 to 0.95 (2) without causing any problem. This restrictive characteristic makes the medium especially suitable for foods. Peptone provides carbon, nitrogen, vitamins and minerals. Dextrose (Glucose) is a carbohydrate source. Phosphate buffers the medium. Magnesium sulfate provides divalent cations and sulfate. Dichloran is an antifungal agent, added to the medium to reduce colony diameters of spreading fungi. This medium can also be used for isolation of fungi from clinical samples. Chloramphenicol is included to inhibit the growth of bacteria present in environmental and food samples. Inhibition of growth of bacteria and restriction of spreading of more-rapidly growing moulds aids in the isolation of slow-growing fungi by preventing their overgrowth by more-rapidly growing species.

Process 40 gm of food sample in a stomacher by adding 200 ml of 0.1% Peptone Water (M028). Shake periodically for 30 minutes with 0.1% Peptone Water for powdered products. Dilute the sample to 1:10 in 0.1% Peptone water and spread on plate. Count the number of Xerophilic colonies per gram of food. The medium can also be used as general medium for the isolation of yeasts and moulds from foodstuffs (2).

Type of specimen

Clinical samples - skin scrapings/ skin lesions, etc.; Food samples

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precaution

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.

Please refer disclaimer Overleaf.

2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further serological and biochemical testing is required for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow granulated free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Medium amber coloured, clear to slightly opalescent gel forms in Petri plates

Reaction

Reaction of 3.16% w/v aqueous solution 22 grams of glycerol at 25°C. pH : 5.6±0.2

pH

5.40-5.80

Cultural Response

Cultural characteristics observed with added 22 grams of glycerol after an incubation at 25°C for upto 6 days.

Organism	Inoculum (CFU)	Growth	Recovery
** <i>Bacillus spizizenii</i> ATCC 6633 (00003*)	≥10 ⁴	inhibited	0%
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	good-luxuriant	≥50%
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 ⁴	inhibited	0%
<i>Mucor racemosus</i> ATCC 42647	50-100	good-luxuriant	
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	50-100	good-luxuriant	≥50%

Key : (*) Corresponding WDCM numbers.

Formerly known as *Bacillus subtilis* subsp. *spizizenii*

Storage and Shelf Life

Store between 15-25°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,4).

Reference


1. Hocking A.D. and Pitt J.I., 1980, J. Appl. Environ. Microbiol., 39:488.
2. Beckers H.J., et al, 1982, Intern. Stand. Org. Document ISO/TC34/SC9/N151.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
5. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC,WagleIndustrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

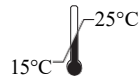

Storage temperature

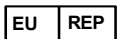

AR Experts BV
Boeingavenue 209
1119 PD Schiphol-Rijk
The Netherlands



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.