



## **Kimmig Fungi Agar Base, Granulated<sup>®</sup>**

**GM1232**

### **Intended Use:**

Recommended for cultivation, isolation and identification of fungi.

### **Composition\*\***

<b>Ingredients</b>	<b>g / L</b>
Peptone	15.000
Sodium chloride	1.000
Dextrose (Glucose)	19.000
Cycloheximide	0.400
Agar	15.000
Final pH ( at 25°C)	6.5±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### **Directions**

Suspend 50.40 grams in 1000 ml purified/distilled water, containing 5ml glycerol. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15lbs pressure (121°C) for 15 minutes. Cool to 45-50°C and aseptically add reconstituted contents of two vials of Kimmig Selective Supplement (FD111) or two vials of George Kimmig Selective Supplement (FD112). Mix well and pour into sterile Petri plates.

### **Principle And Interpretation**

Kimmig Fungi Agar is prepared as described by Kimmig and Rieth (1) for cultivation, isolation, identification and strain preservation of fungi. Fungi identification is usually carried out by examining the hyphae or spores formed by fungi on the medium plates. Rieth later stated that this medium promotes the development of growth forms, which are used as important characteristic criteria in identification (2).

The medium contains peptone, which provides the necessary nitrogenous and carbonaceous nutrients, long chain amino acids, vitamins for the growth of fungi. Dextrose is the fermentable carbohydrate and energy source. Glycerol serves as the carbon source. Kimmig Fungi Agar Base is used as a base for preparation of selective agars for isolation of fungi from heavily contaminated materials. George et al (3) suggested addition of cycloheximide, penicillin and streptomycin while Hantschke (4) suggested the use of colistin and novobiocin.

### **Type of specimen**

Clinical samples - skin or nail scraping or infected hair etc

### **Specimen Collection and Handling:**

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### **Warning and Precautions**

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### **Limitations :**

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement. Further biochemical and serological tests must be carried out for further identification.

### **Performance and Evaluation**

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Cream to yellow colored granular medium.

### Gelling

Firm, comparable with 1.5% Agar gel

### Colour and Clarity of prepared medium

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates

### Reaction

Reaction of 5.04% w/v aqueous solution at 25°C. pH : 6.5±0.2

### pH

6.30-6.70

### Cultural Response

Cultural characteristics observed with added Kimmig Supplement (FD111) or George Kimmig Selective Supplement (FD112), after an incubation at 25- 30°C for 48-72 hours.

### Organism

### Growth

#*Aspergillus brasiliensis* luxuriant

ATCC 16404 (00053\*)

*Candida albicans* ATCC luxuriant

10231 (00054\*)

*Pencillium notatum* ATCC luxuriant

10108

##*Trichophyton* luxuriant

*interdigitale* ATCC 9533

Key: \*Corresponding WDCM numbers.

## Formerly known as *Trichophyton mentagrophytes*

# Formerly known as *Aspergillus niger*

## Storage and Shelf Life

Store dehydrated and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

## Reference

1. Kimmig J. and Rieth H., 1953, *Arzneimittelforsch*, 3:267.
2. Rieth H., 1969, *Mykosen*, 12: 73.
3. George L. K., Ajello L. and Papageorge C., 1954, *J. Lab. Clin. Med.*, 44:422.
4. Hantschke D., 1968, *Mykosen*, 11:769.
5. Isenberg, H.D. *Clinical Microbiology Procedures Handbook* 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) *Manual of Clinical Microbiology*, 11th Edition. Vol. 1.


Revision :08/2026



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



**In vitro diagnostic  
medical device**

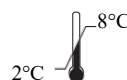

**Storage temperature**



AR Experts BV  
Boeingavenue 209  
1119 PD Schiphol-Rijk  
The Netherlands



**CE Marking**



**Do not use if  
package is damaged**

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.