

## Nutrient Agar No.2, Granulated<sup>®</sup>

GM1269

### Intended Use:

Recommended as a general purpose culture medium, also can be used for microbiological analysis of water as per Czech standard.

### Composition\*\*

Ingredients	g / L
Peptone	10.000
HM peptone B #	10.000
Sodium chloride	5.000
Agar	15.000
Final pH ( at 25°C)	7.2±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# - Equivalent to Beef extract

### Directions

Suspend 40.0 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and aseptically pour into sterile Petri plates.

### Principle And Interpretation

Nutrient Media are general purpose media used for the examination of water and dairy products according to Standard Methods for the Examination of Water and Waste water (1) and Dairy Products (2). Nutrient Agar No. 2 can be used for the microbiological analysis of water as per Czech Standards. It can also be used for sterility testing of aerobes and also for maintenance of subcultures (3). As a general purpose can be used for clinical samples also.

HM peptone B and peptone provide the necessary nitrogen compounds, carbon, long chain amino acids, vitamins and also some trace ingredients for the growth of bacteria. Sodium chloride maintains osmotic equilibrium of the medium.

### Type of specimen

Clinical samples- faeces, urine, etc.; Food and dairy samples; Water samples

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5). For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (2,6,7). For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (1). After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations:

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further serological and biochemical testing is required for complete identification.

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Cream to yellow homogeneous granular media.

### Gelling

Firm, comparable with 1.5% Agar gel

### Colour and Clarity of prepared medium

Light yellow to amber coloured clear to slightly opalescent gel forms in Petri plates

### Reaction

Reaction of 4.0% w/v aqueous solution at 25°C. pH : 7.2±0.2

### pH

7.00-7.40

### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours

Organism	Inoculum (CFU)	Growth	Recovery
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	≥70%
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	luxuriant	≥70%
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	50-100	luxuriant	≥70%
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	luxuriant	≥70%

Key : (\*) Corresponding WDCM numbers (#) Formerly known as *Enterobacter aerogenes*.

## Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

## Reference

- Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
- American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
- MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Lippincott, Williams and Wilkins, Baltimore.
- Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
- Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
- Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 06/2026




HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC,WagleIndustrial Area,  
Thane (W) -400604, MS, India



**IVD** *In vitro* diagnostic  
medical device

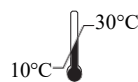

**Storage temperature**

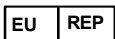

AR Experts BV  
Boeingavenue 209  
1119 PD Schiphol-Rijk  
The Netherlands



**CE Marking**



**Do not use if  
package is damaged**

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory,diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.