



Malt Extract Agar Base (w/ Mycological Peptone), Granulated[®]

GM137

Intended Use:

For detection, isolation and enumeration of yeasts and moulds.

Composition**

Ingredients	g / L
Malt extract	30.000
Mycological peptone	5.000
Agar	15.000
Final pH (at 25°C)	5.4±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 50.0 grams in 1000 ml purified/distilled water and soak for 15 minutes. Sterilize by autoclaving at 115°C(10 lbs pressure) for 10 minutes. Mix well before dispensing. Avoid overheating. If desired, to adjust acidic pH use 10% Lactic Acid (FD095).

Principle And Interpretation

The laboratory diagnosis of fungal infection relies largely on direct as opposed to indirect methods. The use of malt and malt extracts for the propagation of yeasts and moulds is quite common. Reddish (1) described a culture medium prepared from malt extract that was a satisfactory substitute for wort. Malt Extract Medium is similar to the formula of Galloway and Burgess (2) used for the detection, isolation and enumeration of yeasts and moulds.

Malt extract provides an acidic environment and nutrients favorable for growth and metabolism of yeasts and moulds. Mycological peptone rapidly gives a luxuriant growth with typical morphology and pigmentation. For mycological count, it is advisable to adjust the reaction of medium more acidic with addition of 10% lactic acid. Antibiotics may be added as sterile solutions to the molten medium immediately before pouring into sterile Petri plates (3) in order to suppress bacterial growth.

Type of specimen

Clinical samples : skin scrapings. nail scrapings etc.

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further biochemical tests must be carried out for further identification

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to beige coloured granulated media

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Amber coloured clear to slightly opalescent gel forms in Petri plates

Reaction

Reaction of 5.0% w/v aqueous solution at 25°C. pH : 5.4±0.2

pH

5.20-5.60

Cultural Response

Cultural characteristics observed after an incubation at 25-30°C for 48-72 hours. Incubate for 7 days for Trichophyton species.

Organism	Inoculum (CFU)	Growth	Recovery
# <i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	50-100	luxuriant	-
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	luxuriant	≥70%
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	50-100	luxuriant	≥70%
<i>Penicillium chrysogenum</i> ATCC 9179	50-100	luxuriant	-
## <i>Trichophyton interdigitale</i> ATCC 9533	50-100	luxuriant	-

Key : *Corresponding WDCM numbers.

Formerly known as *Aspergillus niger*

Formerly known as *Trichophyton mentagrophytes*

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference

1. Reddish A., 1919, Abstr. Bacteriol., 3:6.
2. FDA Bacteriological Analytical Manual, 2005, 18th Ed., AOAC, Washington, DC.
3. Galloway L. D. and Burgess R., 1952, Applied Mycology and Bacteriology, 3rd Ed., Leonard Hill, London, pg. 54 and 57.
4. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 07/2026



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC,Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

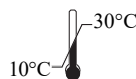

Storage temperature



AR Experts BV
Boeingavenue 209
1119 PD Schiphol-Rijk
The Netherlands



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.