



Tryptose Broth, Granulated[®]

GM177

Intended Use:

Recommended for the cultivation primarily of *Brucella* species.

Composition**

Ingredients	g / L
Tryptose	20.000
Glucose(Dextrose)	1.000
Sodium chloride	5.000
Final pH (at 25°C)	7.3±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 26.0 grams in 1000 ml purified/distilled water. If desired, add 0.5-1.0% agar to the medium. Heat if necessary to dissolve the medium completely. Dispense in tubes or flasks as desired. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle And Interpretation

Huddleson used Tryptose media for the isolation of *Brucella* species from man (1). Tryptose containing media, rather than the conventionally used meat infusion media have been used for the enumeration and isolation of *Brucella* species (2,3). Tryptose Broth is also recommended by APHA (4) and FDA (5). This medium can be used as general purpose media for cultivation of wide variety of organisms. It can also be supplemented with defibrinated blood (sheep, horse) to prepare blood containing medium for the isolation of fastidious organisms like *Brucella*. Tryptose Broth can be supplemented with 0.1% agar for the cultivation of anaerobes. Dextrose is the source of energy. Tryptose serves as nitrogen source while sodium chloride maintains osmotic equilibrium.

Type of specimen

Clinical sample : CSF; Food samples

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,7).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (4,5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. All presumptive anaerobic organisms must be identified by confirmatory test.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow colored granular medium

Colour and Clarity of prepared medium

Basal Medium : Yellow coloured, clear solution. With addition of 5% v/v sterile defibrinated blood, cherry red coloured, opaque solution forms in tubes.

Reaction

Reaction of 2.6% w/v aqueous solution at 25°C. pH : 7.3±0.2

pH

7.10-7.50

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 48-72 hours with added 5% v/v sterile defibrinated blood in presence of 10% Carbon dioxide (CO₂).

Organism	Inoculum (CFU)	Growth
<i>Brucella melitensis</i> ATCC 4309	50-100	good-luxuriant
<i>Brucella suis</i> ATCC 4314	50-100	good-luxuriant
<i>Streptococcus pneumoniae</i> ATCC 6303	50-100	good-luxuriant
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant

Storage and Shelf Life

Store below 10-30°C in a tightly closed container and the prepared medium at 15-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,7).

Reference

- Huddleson I. F., 1943, Brucellosis in man and animals, rev., Ed., The Commonwealth Fund, New York, N.Y.
- Huddleson I. F., 1939, Brucellosis in Man and Animals, Oxford University Press, Oxford, England.
- Ruiz Castañeda M., 1947, Proc. Soc. Exp. Biol. Med., 64:114.
- Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- U.S. Food and Drug Administration, 1995, Bacteriological Analytical Manual, 8th Ed., AOAC International, Gaithersburg, Md.
- Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
- Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W.(2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 03/2026



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

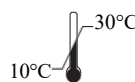

Storage temperature



AR Experts BV
Boeingavenue 209
1119 PD Schiphol-Rijk
The Netherlands



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.