

## Lactobacillus MRS Broth, Granulated (MRS Broth, Granulated) GM369

### Intended Use

Recommended for cultivation of all Lactobacilli from clinical and non clinical samples.

### Composition\*\*

Ingredients	g / L
Proteose peptone	10.000
HM Peptone B#	10.000
Yeast extract	5.000
Dextrose(Glucose)	20.000
Polysorbate 80 (Tween 80)	1.000
Ammonium citrate	2.000
Sodium acetate	5.000
Magnesium sulphate	0.100
Manganese sulphate	0.050
Dipotassium hydrogen phosphate	2.000
Final pH ( at 25°C)	6.5±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# Equivalent to Beef Extract

### Directions

Suspend 55.15 grams in 1000 ml purified/distilled water. Heat if necessary to dissolve the medium completely. Distribute in tubes, bottles or flasks as desired. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

### Principle And Interpretation

Lactobacilli MRS media are based on the formulation of deMan, Rogosa and Sharpe (1) with slight modification. It supports luxuriant growth of all Lactobacilli from oral cavity (1), dairy products (2), foods (3), faeces (4,5) and other sources (6).

Proteose peptone and HM peptone B supply nitrogenous and carbonaceous compounds. Yeast extract provides vitamin B complex and dextrose is the fermentable carbohydrate and energy source. Polysorbate 80 supplies fatty acids required for the metabolism of Lactobacilli. Sodium acetate and ammonium citrate inhibit Streptococci, moulds and many other microorganisms.

### Type of specimen

Clinical samples - faeces, swab from oral cavity; Food and dairy samples

### Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5). For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (2,7,8). After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations :

1. Due to nutritional variation, some strains may show poor growth.
2. Further biochemical and serological tests must be carried out for complete identification.

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Cream to yellow colored granular medium

### Colour and Clarity of prepared medium

Light amber coloured, clear to slightly opalescent solution in tubes

### Reaction

Reaction of 5.51% w/v aqueous solution at 25°C. pH : 6.5±0.2

### pH

6.30-6.70

### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours or longer (with 5% CO<sub>2</sub>)

Organism	Inoculum (CFU)	Growth
<i>Lactobacillus fermentum</i> ATCC 9338	50-100	luxuriant
<i>Lactobacillus leichmannii</i> ATCC 7830	50-100	luxuriant
<sup>^</sup> <i>Lactiplantibacillus</i> <i>plantarum</i> ATCC 8014	50-100	luxuriant
<i>Lactobacillus rhamnosus</i> ATCC 9595	50-100	luxuriant
<i>Lactobacillus sakei</i> ATCC 15521 (00015*)	50-100	luxuriant
<i>Lactobacillus lactis</i> ATCC 19435 (00016*)	50-100	luxuriant
<i>Pediococcus pentosaceus</i> ATCC 33316 (00158*)	50-100	luxuriant

Key: (\*) Corresponding WDCM numbers.

<sup>^</sup> Formerly known as *Lactobacillus plantarum*

## Storage and Shelf Life

Store dehydrated and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

## Reference


1. deMan J., Rogosa M. and Sharpe M., 1960, J. Appl. Bacteriol., 23:130.
2. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
3. Marshall R.T. (Ed.), 1992, Standard Methods for the Examination of Dairy Products, 16th ed., APHA, Washington, D.C.
4. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
6. MacFaddin J., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol.1, Williams and Wilkins, Baltimore.
7. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
8. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



**In vitro diagnostic  
medical device**

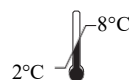

**Storage temperature**



CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
[www.cepartner4u.eu](http://www.cepartner4u.eu)



**CE Marking**



**Do not use if  
package is damaged**

#### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.