

## Magnetic Tip Comb for Insta NX<sup>®</sup> Mag 96

<u>Product Name</u>	<u>Product Code</u>	<u>Packing</u>
Magnetic Tip Comb for Insta NX <sup>®</sup> Mag96	LA1097A	40 NO



Fig.1 Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96

### Intended Use

Recommended for adsorption and release of magnetic beads in the sample for purification of nucleic acids.

**Introduction:** The Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 is used for purification of nucleic acids. It helps magnetic beads to adsorb, transfer and release in the sample. The Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 is used for a processing 96 samples at a time. The Magnetic Tip Comb is compatible for use with solutions containing ethanol or isopropanol.


### Specifications:

The Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 is compatible with Insta NX<sup>®</sup> Mag96 (LA1097).

**Storage conditions:** The Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 can be stored at Room temperature (15-25°C).

### Directions for use:

1. Slide open the door of the Insta NX<sup>®</sup> Mag96.
2. Select the desired location on the screen of Insta NX<sup>®</sup> Mag96 for placing Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 in the machine.
3. Place the Magnetic Tip Comb for Insta NX<sup>®</sup> Mag96 at the desired position on the machine.
4. Place the filled 96 Deep Well Plates onto the machine.
5. Close the door of the Insta NX<sup>®</sup> Mag96.

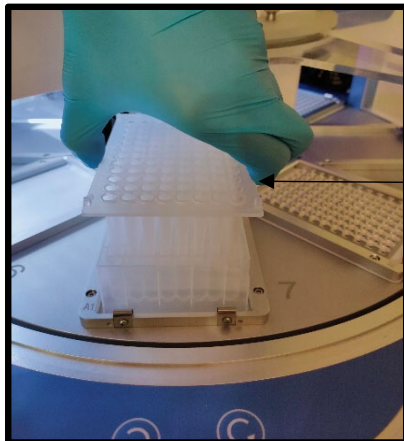

**Magnetic Tip Comb for Insta NX®  
Mag96**

Fig.2 Inserting Magnetic Tip Comb in Insta NX® Mag96

### **Disposal**

User must ensure safe disposal by autoclaving and/ or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed off in accordance with current laboratory techniques.

### **Warning and Precautions**

Not for Medicinal Use. Read the procedure carefully. Wear protective gloves/protective clothing/eye protection/face protection. Follow good clinical laboratory practices while handling clinical samples. Standard precautions should be followed as per established guidelines. Safety guidelines may be referred in safety data sheets of the product.

### **Safety Information**


Take appropriate laboratory safety measures and wear gloves when handling. Not compatible with disinfecting agents containing bleach. Please refer the Safety Data Sheet (SDS) for information regarding hazards and safe handling practices.

### **Technical Assistance**

At HiMedia, we pride ourselves on the quality and availability of our technical support. For any kind of technical assistance, mail to [mb@himedialabs.com](mailto:mb@himedialabs.com).

---

Please refer disclaimer Overleaf.

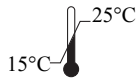

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Private Limited,  
Reg. Off: Plot No. C-40, Road No. 21Y,  
MIDC, Wagle Industrial Area, Thane,  
(West) 400604, Maharashtra, INDIA.  
Web: [www.himedialabs.com](http://www.himedialabs.com)



07/2025

PILA1097A\_0/0722

LA1097A-00

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. Reg.office : Plot No. C-40, Road No. 21Y, MIDC, Wagle Industrial Area, Thane, (West) 400604, Maharashtra, INDIA.  
Customer Care No.: 00-91-22-6116 9797 Tel: 00-91-22-6147 1919, 6903 4800 Email: [techhelp@himedialabs.com](mailto:techhelp@himedialabs.com) Website: [www.himedialabs.com](http://www.himedialabs.com)