



## Sabouraud Maltose Agar

M062

### Intended Use:

Recommended for propagation of yeasts and moulds, particularly the parasitic fungi concerned with skin and scalp lesions.

### Composition\*\*

Ingredients	g / L
Maltose	40.000
Mycological peptone	10.000
Agar	15.000
Final pH ( at 25°C)	5.6±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 65 grams in 1000 ml purified / distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

Fungi were among the first microorganisms recognized because some of the fruiting structures, such as the mushrooms, are large enough to be seen without a microscope. Fungi can be grouped simply on the basis of morphology as either yeasts or moulds (1). Sabouraud Maltose Agar was formulated by Sabouraud (2) and is used for the isolation and differentiation of yeast and moulds. (3,4,5)

Mycological peptone provides nitrogen, vitamins, minerals, amino acids and growth factors. Maltose provides an energy source for the growth of microorganisms. The low pH favours fungal growth and inhibits contaminating bacteria from clinical specimens (6). The acid reaction of the final medium is inhibitory to a large number of bacteria making it particularly useful for cultivating fungi and aciduric microorganisms. For isolation of fungi from contaminated specimens, a selective medium should be inoculated simultaneously. Incubate cultures for 4 to 6 weeks before reporting as negative.

### Type of specimen

Clinical samples - skin and scalp lesions, etc.

### Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,7).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement. Further biochemical and serological tests must be carried out for further identification.

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Cream to yellow homogeneous free flowing powder

**Gelling**

Firm, comparable with 1.5% Agar gel

**Colour and Clarity of prepared medium**

Light amber coloured clear to slightly opalescent gel forms in Petri plates

**Reaction**

Reaction of 6.5% w/v aqueous solution at 25°C. pH : 5.6±0.2

**pH**

5.40-5.80

**Cultural Response**

Cultural characteristics observed after an incubation at 25 - 30°C for 48 - 72 hours.(Incubate *Trichophyton* species for upto 7 days)

Organism	Inoculum (CFU)	Growth	Recovery
# <i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	50-100	good-luxuriant	
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	good-luxuriant	>=70%
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant (Inhibited on media with lower pH)	>=70%
<i>Lactobacillus casei</i> ATCC 9595	50-100	good-luxuriant	>=70%
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	50-100	good-luxuriant	>=70%
<i>Trichophyton rubrum</i> ATCC 28191	50-100	good-luxuriant	

Key : \*Corresponding WDCM numbers.

#- Formerly known as *Aspergillus niger*

**Storage and Shelf Life**

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

**Disposal**

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,7).

**Reference**

1. Murray P. R., Baron J. H., Pfaller M. A., Tenover J. C. and Tenover F. C., (Ed.), 2003, Manual of Clinical Microbiology, 8th Ed., American Society for Microbiology, Washington, D.C.
2. Sabouraud R., 1892, Ann. Dermatol. Syphil. 3 : 1061.
3. Davidson and Dowding, 1932, Arch. Dermatol. Syphilol. 26:660.
4. Davidson, Dowding and Buller. 1932. Can. J. Res. 6:1.
5. Frank L. S., 1932, Arch. Dermatol. Syphilol., 26: 457
6. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
7. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.

Revision :04/2024



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
[www.cepartner4u.eu](http://www.cepartner4u.eu)



**In vitro diagnostic  
medical device**



**CE Marking**

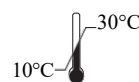

**Storage temperature**



**Do not use if  
package is damaged**

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.