



## KF Streptococcus Agar Base w/ BCP

M1007

### Intended Use

Recommended for detection and enumeration of faecal *Streptococci* from clinical and food samples.

### Composition\*\*

Ingredients	g/ L
Proteose peptone	10.000
Yeast extract	10.000
Sodium chloride	5.000
Sodium $\beta$ -glycerophosphate	10.000
Maltose	20.000
Lactose	1.000
Sodium azide	0.400
Bromocresol purple	0.015
Agar	20.000
Final pH ( at 25°C)	7.2 $\pm$ 0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 76.41 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. **DO NOT AUTOCLAVE.** Overheating will lower the pH and render the medium less productive. Cool to 45-50°C and aseptically add 10 ml of 1% 2,3,5-Triphenyl Tetrazolium Chloride (TTC) (FD057) to sterile medium. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

*Streptococci* are spherical, gram-positive bacteria and form a part of the normal commensal flora of the mouth, skin, intestine, upper respiratory tract of humans. Streptococci found in the faeces form the faecal *Streptococci* and constitute of *Streptococci* with group D Lancefield antigens. The types include *Streptococcus faecalis*, *Streptococcus faecium*, *Streptococcus bovis* and *Streptococcus duran*. They are low-grade pathogens and rarely cause disease. However, they may cause urinary tract infection in catheterized patients; mixed abdominal wound infections following gut surgery; and endocarditis on abnormal valves. Kenner - Faecal (KF) Medium were developed by Kenner et al (1,2) for detecting Streptococci in water and food materials.

Proteose peptone along with yeast extract provide nitrogen, carbon, sulphur, amino acids, vitamins and trace ingredients to the faecal Streptococci. Lactose and maltose are the fermentable carbohydrates and therefore serve as energy sources. Sodium azide is a selective agent, which hampers the growth of gram-negative bacteria.

2,3,5-Triphenyl Tetrazolium Chloride is reduced to insoluble formazan by actively metabolizing cells, resulting in the formation of pink or red colonies. Bacteria resistant to azide, utilize lactose and / or maltose. The acidity so produced changes the colour of the indicator dyes to yellow. Bacterial cells reduce TTC to insoluble formazan, resulting in the formation of pink to red colonies.

### Type of specimen

Clinical samples - faeces, Food and Water samples.

### Specimen Collection and Handling:

Samples can be directly streaked or sterile membrane filters through which the water samples have been passed are aseptically placed on the media. After an incubation at 35-37°C for 24-48 hours, Enterococci appear as pink to red colonies. After this presumptive identification, further confirmatory tests should be carried out (3,4). After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. The media should be handled by trained personnel only. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines

should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations:

1. KF medium was devised to detect fetal contamination in water and is not specific for presumptive identification of the group D *Streptococcus* species, further tests are required along with KF (4).

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Cream to greyish yellow homogeneous free flowing powder

#### Gelling

Firm, comparable with 2.0% agar gel.

#### Colour and Clarity of prepared medium

Light purple coloured, clear to slightly opalescent gel forms in Petri plates

#### Reaction

Reaction of 7.64% w/v aqueous solution at 25°C. pH : 7.2±0.2

#### pH

7.00-7.40

#### Cultural Response

Cultural characteristics observed with added FD057, after an incubation at 35-37°C for 48-72 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	≥10 <sup>4</sup>	inhibited	0%	
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	good-luxuriant	≥50%	red-maroon
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 <sup>4</sup>	inhibited	0%	

Key : (\*) - Corresponding WDCM numbers (#) - Formerly known as *Enterobacter aerogenes*

### Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

### Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

### Reference

1. Kenner B. A., Clark H. F. and Kabler P. W., 1960, Am. J. Public Health, 50:1553.
2. Kenner B. A., Clark H. F. and Kabler P. W., 1961, Appl. Microbiol., 9:15.
3. Facklam R. R. and Moody M. P., 1970, Appl. Microbiol., 20:245.
4. MacFaddin J. F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams and Wilkins, Baltimore.
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
7. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
[www.cepartner4u.eu](http://www.cepartner4u.eu)



**In vitro diagnostic  
medical device**



**CE Marking**

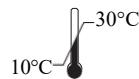

**Storage temperature**



**Do not use if  
package is damaged**

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.