



Candida Medium

M104

Intended Use:

Recommended for selective isolation and cultivation of *Candida* species from clinical and non-clinical samples.

Composition**

Ingredients	g / L
Mycological peptone	2.500
Dextrose (Glucose)	5.000
Disodium hydrogen phosphate	5.000
Sodium sulphite	5.000
Bismuth sulphite indicator	3.000
Agar	15.000
Final pH (at 25°C)	7.6±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 35.5 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. **DO NOT AUTOCLAVE**. Cool to 45-50°C and aseptically add 0.3 units of Penicillin and 25 mcg Streptomycin per ml of sterile medium. Mix well and pour into sterile Petri plates.

Principle And Interpretation

Candida is a genus of yeasts responsible for infections such as oropharyngeal candidiasis, vaginal candidiasis and candidemia. Candida Medium is used for the selective cultivation and differentiation of *Candida* species. Candida Medium was originally developed by Nickerson (1). It is also used for processing and inoculation of specimens like tissues, skin scraping, nails and hair (2,3).

Mycological peptone in the medium provides essential nitrogenous nutrients while dextrose acts as carbon source and phosphate maintains buffering action of medium. This medium also contains sodium sulphite, which is reduced by *Candida* species to form sulphide. Bismuth in the medium combines with the sulphide to produce brown to black pigmented colonies and zones of dark precipitate in the medium surrounding the colonies of some species. Bismuth sulphite also acts as an inhibitor of bacterial growth. Selectivity of medium is increased by incorporation of penicillin and streptomycin in the medium, which helps to suppress the growth of many bacteria.

Differentiation of *Candida* is based on the growth patterns and pigmentation of isolated colonies.

Type of specimen

Clinical samples - skin scraping, nails and hair

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1.This medium is general purpose medium and may not support the growth of fastidious organisms.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Light yellow coloured, clear to slightly opalescent gel forms in Petri plates

Reaction

Reaction of 3.55% w/v aqueous solution at 25°C. pH : 7.6±0.2

pH

7.40-7.80

Cultural Response

Cultural characteristics observed after an incubation at 30°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	good-luxuriant	≥50%
<i>Candida tropicalis</i> ATCC 1369	50-100	good-luxuriant	≥50%
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 ⁴	inhibited	0%

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference

1. Nickerson, 1953, J. Infect. Dis., 93:43.
2. Emmons, Binford, Utz and Kwon-Chung, 1977, Medical Mycology, 3rd Ed. , W. B. Saunders Co., Philadelphia.
3. Haley, Trandel and Coyle, 1980, Cumitech 11, Practical Methods for Culture and Identification Of Fungi In The Clinical Mycology Laboratory, Coord Ed., Sherris, ASM, Washington, D.C.
4. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision :04/ 2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

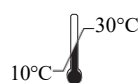

Storage temperature

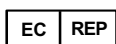

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.