



B.T.B. Lactose Agar, Modified (Lactose Blue Agar)

M1081

Intended Use:

Differentiation of lactose fermenting and non-fermenting bacteria belonging to *Enterobacteriaceae*.

Composition**

Ingredients	g / L
Peptone	3.500
Tryptone	3.500
Sodium chloride	5.000
Lactose	15.500
Bromo thymol blue	0.040
Agar	13.000
Final pH (at 25°C)	7.0±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 40.54 grams in 1000 ml purified / distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

Principle And Interpretation

Reactions with lactose are of great practical importance for the primary isolation of *Enterobacteria* from clinical specimens. The specimens e.g. faeces is usually plated on a lactose-containing medium on which lactose fermenters and lactose nonfermenters form coloured and pale colonies respectively due to the dye incorporated. This procedure makes an immediate presumptive distinction between colonies of the true intestinal pathogens possible. *Salmonella* and *Shigella*, do not ferment lactose while the common intestinal commensals, *Escherichia* and *Klebsiella*, which do ferment lactose (1). Lactose Blue Agar is used for differentiating lactose fermenting and non-fermenting bacteria belonging to the family *Enterobacteriaceae*.

Tryptone and peptone provide essential nutrients for bacterial metabolism. Lactose provides a fermentable carbohydrate source for the enteric bacteria. Bromo thymol blue is the pH indicator for indicating acid production due to carbohydrate fermentation. The dye turns yellow at acidic pH and imparts yellow colour to the colony. Alkalinization produces a blue coloration. Winkle (2) recommended addition of 0.28g/l metachrome yellow to suppress the swarming of *Proteus* species.

Type of specimen

Clinical samples: Isolated Microorganisms from clinical samples, faeces

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further biochemical testing is required for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to greenish yellow homogeneous free flowing powder

Gelling

Firm, comparable with 1.3% Agar gel.

Colour and Clarity of prepared medium

Green coloured, clear to slightly opalescent gel forms in Petri plates.

Reaction

Reaction of 4.05% w/v aqueous solution at 25°C. pH : 7.0±0.2

pH

6.80-7.20

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of Colony
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	≥70%	yellow, opaque
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	luxuriant	≥70%	bluish
<i>Salmonella</i> Typhi ATCC 6539	50-100	luxuriant	≥70%	bluish
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	good-luxuriant	≥70%	deep yellow

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,4).

Reference

1. Cruikshank R., Duguid J. P., Marmion B. P., Swain R. H. A., (Eds.), 1975, Medical Microbiology, The Practice of Medical Microbiology, 12th Edition, Vol. II, Churchill Livingstone.
2. Winkle S., 1947, Zbl. Bakt. I. Orig., 152:103.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.

Revision : 06/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



In vitro diagnostic
medical device



CE Marking

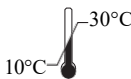

Storage temperature



Do not use if
package is damaged

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.