



## H.S. Vaccine Medium (Standard Nutrient Broth)

M116

### Intended Use:

Used for the large scale cultivation of bacteria for production of vaccine preparations.

### Composition\*\*

Ingredients	g / L
Peptone	10.000
HM infusion, solid#	10.000
Sodium chloride	5.000
Final pH ( at 25°C)	7.6±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# Equivalent to Peptic digest of lean meat infusion, solids

### Directions

Suspend 25 grams in 1000 ml purified / distilled water. Heat if necessary to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C.

### Principle And Interpretation

Peptone is the principal source of organic nitrogen while HM infusion, solid provides carbohydrates, vitamins, organic nitrogen compounds and salts. Sodium chloride maintains osmolality of the medium.

### Type of specimen

Isolated Microorganism from clinical samples.

### Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (1,2).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations :

1.This medium is general purpose medium and may not support the growth of fastidious organisms.

### Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Cream to light yellow homogeneous free flowing powder

#### Colour and Clarity of prepared medium

Amber coloured clear solution without any precipitate.

#### Reaction

Reaction of 2.5% w/v aqueous solution at 25°C. pH : 7.6±0.2

#### pH

7.40-7.80

#### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-48 hours.

Organism	Inoculum (CFU)	Growth
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	good-luxuriant
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	good-luxuriant
<i>Staphylococcus epidermidis</i> ATCC 12228 (00036*)	50-100	good-luxuriant
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant

Key : (#) *Enterobacter aerogenes* ATCC 13048

(\*) Corresponding WDCM numbers

### Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

### Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (1,2).

### Reference

1. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2<sup>nd</sup> Edition.
2. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 04/2024




HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



**In vitro diagnostic  
medical device**

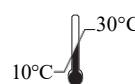

**Storage temperature**

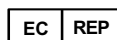

CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
www.cepartner4u.eu



**CE Marking**



**Do not use if  
package is damaged**

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.