



Yersinia Selective Broth Base

M1861

Intended use

Recommended for the selective enrichment of *Yersinia enterocolitica*.

Composition**

Ingredients	g / L
Peptone	10.000
L-asparaginic acid	20.000
Sodium pyruvate	2.500
Tween 80	2.200
MOPS/TRIS	5.500
Final pH (at 25°C)	7.2±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 40.2 grams in 1000 ml purified / distilled water. Heat if necessary to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C and aseptically add the rehydrated contents of 1 vial of BT Selective Supplement (FD286). Mix well and dispense into sterile tubes or flasks as desired.

Principle And Interpretation

Yersinia enterocolitica is widely distributed in lakes and reservoirs. Epizootic outbreaks of diarrhea, lymphadenopathy, pneumonia and spontaneous abortions occur in various animals. It is the most common species of *Yersinia* recovered from clinical specimens. *Yersinia enterocolitica* is biochemically more active at room temperature than at 37°C. *Yersinia* Selective Broth Base is recommended for the selective enrichment of *Yersinia enterocolitica* (1).

The medium contains peptone which provides carbonaceous and nitrogenous compounds, long chain amino acids and other essential compounds. L-asparaginic acid and sodium pyruvate are growth supporting ingredients. MOPS/TRIS buffers the medium.

Type of specimen

Clinical samples - stool sample; Food and dairy samples

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (2,3).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (4,5,6).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Further serological and Biochemical tests must be carried out on pure colony for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow coloured, homogeneous free flowing powder

Colour and Clarity of prepared medium

Yellow coloured clear solution without any precipitate.

Reaction

Reaction of 4.02 % w/v aqueous solution at 25°C . pH : 7.2±0.2

pH

7.00-7.40

Cultural Response

Cultural characteristics observed after an incubation at 25-30°C for 18- 24 hours.

Organism	Inoculum (CFU)	Growth
<i>Yersinia enterocolitica</i> ATCC 27729	50-100	good-luxuriant
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	≥10 ⁴	inhibited
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	≥10 ⁴	inhibited

Key :(*)- Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (2,3).

Reference

1. Schiemann D. A., 1979, Can. J. Microbiol., 25: 1298.
2. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
3. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
4. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
5. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 06/2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

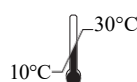

Storage temperature

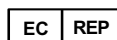

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.