



Technical Data

Buffered Peptone Water (6 Fold Strength phosphate buffer)

M2050

Intended Use:

Recommended as a pre-enrichment medium used for increasing the recovery of injured *Salmonella* species from foods prior to selective enrichment and isolation. It can also be used for clinical samples.

Composition**

Ingredients	Gms / Litre
Tryptone	10.000
Sodium chloride	5.000
Disodium hydrogen phosphate, 12H ₂ O	54.000
Potassium dihydrogen phosphate	9.000
Final pH (at 25°C)	7.0±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 45.41 grams (equivalent weight of dehydrated medium per litre) in 1000 ml purified / distilled water. Heat if necessary to dissolve the medium completely. Dispense in tubes or flasks as desired. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes.

Principle And Interpretation

Buffered Peptone Water (6 fold strength phosphate buffer) is a pre-enrichment medium designed to help recovery of sub-lethally damaged *Salmonellae* before transfer to a selective medium. This pre-enrichment medium is free from inhibitors and is well buffered and provides conditions for resuscitation of the cells that have been injured by processes of food preservation.

The media contain tryptone as a source of carbon, nitrogen, long chain amino acids, vitamins, minerals and other essential nutrients. Sodium chloride maintains the osmotic balance and phosphates buffer the medium. The broth is rich in nutrients and produces high resuscitation rates for sub-lethally injured bacteria and supports intense growth. The phosphate buffer system prevents bacterial damage due to changes in the pH of the medium.

Type of specimen

Clinical samples: Pathological specimens, Food and dairy samples, Water samples.

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (1,5,6).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (2).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Though this medium is selective for *Salmonella* other species of *Enterobacteriaceae* may grow.

2. *Salmonella* Typhi and *Shigella* species may not grow on this medium.
3. Moreover *Proteus*, *Pseudomonas* and *Citrobacter* species may mimic enteric pathogens by producing small red colonies.
4. Final confirmation of suspected colonies must be carried out by serological and biochemical tests.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Colour and Clarity of prepared medium

Cream to light yellow coloured, clear solution without any precipitate

Reaction

Reaction of 4.54% w/v aqueous solution at 25°C. pH : 7.0±0.2

pH

6.80-7.20

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.(Recovery is carried out using XLD Agar, M031).

Organism	Inoculum (CFU)	Growth	Recovery
<i>Salmonella</i> Enteritidis ATCC 50-100 13076 (00030*)		good-luxuriant	≥50%
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	≥50%
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	good-luxuriant	≥50%

Key: (*) Corresponding WDCM numbers

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 15-25°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,4).

Reference

1. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
2. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
5. Salfinger Y., and Tortorello M.L. 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 01 / 2019



In vitro diagnostic medical device



CE Marking

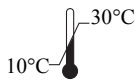

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U, Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.