



McBride Listeria Agar Base

M386

Intended Use:

Recommended for selective isolation and cultivation of *Listeria monocytogenes* from foodstuffs, clinical samples.

Composition**

Ingredients	g / L
Tryptose	10.000
HM peptone B #	3.000
Sodium chloride	5.000
Glycine anhydride	10.000
Lithium chloride	0.500
Phenyl ethanol	2.500
Agar	15.000
Final pH (at 25°C)	7.3±0.2

**Formula adjusted, standardized to suit performance parameters

Equivalent to Beef extract

Directions

Suspend 46.0 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Before gelling, aseptically add sterile defibrinated blood to a final concentration of 5% v/v and add rehydrated contents of one vial of McBride Selective Supplement (FD070). Mix well and pour into sterile Petri plates.

Principle And Interpretation

The disease listeriosis is a frequent cause of abortions in cattle and sheep. In human, symptoms are manifested as septicemia, encephalitis and circulatory monocytosis (1). *Listeria* multiplies over a wide range of temperatures, from 3°C to 45°C, and over a pH range of 5.0 to 9.6. It also survives in food products with pH levels outside these parameters (2). Because of these properties, *Listeria* survives the various food processing techniques (3). McBride Listeria Agar (4), recommended by APHA (5) is used for isolating *Listeria* from clinical specimens and foodstuffs including raw milk (6). This medium helps in the detection of low numbers of *L. monocytogenes* present in food samples. Tryptose and HM peptone B in the medium supply nitrogen, carbon, sulphur and trace nutrients required for the growth of *Listeria*. Phenyl ethyl alcohol is bacteriostatic for gram-negative bacteria as it selectively inhibits DNA synthesis (7). Sodium chloride maintains the osmotic balance of the medium. Glycine inhibits certain gram-negative and gram-positive bacteria including *Escherichia coli* and *Enterococcus faecalis*, the common accompanying contaminants. Lithium chloride also has antibacterial activity. Further selectivity is achieved by the addition of McBride Selective Supplement (FD070). The detection of *L. monocytogenes* is greatly improved by pre-enrichment in liquid media either by one step or two steps. In one step method (8), infected material is inoculated directly in Listeria Selective Broth Base (M889), while in two steps method (9), infected material is inoculated in Listeria Enrichment Broth Base (UVM) (M890A) and incubated at refrigeration temperature of 4°C for few weeks (cold enrichment), as the organism has the ability to grow in low temperature. It is then inoculated in Fraser Secondary Enrichment Broth Base (M1083), followed by plating onto selective agar such as McBride Listeria Agar. The presumptive *Listeria* colonies are selected under 45° transillumination. *Listeria* colonies are dense white to iridescent white appearing as crushed glass. Small colonies tend to be blue, while non-*Listeria* show yellowish orange colonies that are further purified. McBride Listeria Agar can be used as a plating medium with or without supplementation of blood.

Type of specimen

Clinical samples - Body tissues or body fluids, vaginal swabs ; Dairy samples

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (10,11).

For dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (6).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. Further biochemical tests are needed for a final identification of the isolated organisms.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Basal medium : Light amber coloured clear to very slightly opalescent gel. After addition of 5%v/v sterile blood : Cherry red opaque gel forms in Petri plates

Reaction

Reaction of 4.6% w/v aqueous solution at 25°C. pH : 7.3±0.2

pH

7.10-7.50

Cultural Response

Cultural characteristics observed in anaerobic atmosphere with added McBride Selective Supplement (FD070) and 5%v/v sterile defibrinated blood, after an incubation at 35-37°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth w/ FD070	Recovery w/ FD070	Growth w/ blood and FD070	Recovery w/ blood and FD070
<i>Listeria monocytogenes</i> ATCC 19112	50-100	good-luxuriant	≥50%	good-luxuriant	≥50%
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	none-poor	≤10%	none-poor	≤10%
<i>Pseudomonas aeruginosa</i> ATCC 27853(00025*)	50-100	none-poor	≤10%	none-poor	≤10%
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	none-poor	≤10%	none-poor	≤10%

* Corresponding WDCM Numbers

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (10,11).

Reference

1. Hyslop N., St. G. and Osborne A. D., 1959, Vet. Rec. 71 : 1082
2. Conner D. E., Brackett R. E., and Beuchat L. R., 1986, Appl. Environ. Microbiol., 52: 59
3. Doyle M. P., Meske L. M. and Marth E. H., 1985, J. Food Prot., 48:740
4. McBride M. E. and Girard F., 1960, J. Lab. Clin. Med., 55:153.
5. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, American Public Health Association, Washington, D.C.
6. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
7. Dowell, Hill and Altemeier, 1964, J. Bacteriol., 88:1811.
8. FDA Bacteriological Analytical Manual, 2005, 18th Ed., AOAC, Washington, DC.
9. McClain D. and Lec W., 1989, Lab. Comm No. 57, Revised, May 24, 1989, U.S. Dept. of Agric., FSIS, Microbiol. Div., Beltsville, MD.
10. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
11. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 05/2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

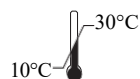

Storage temperature

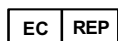

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.