



## Litmus Lactose Bile Salt Agar (LLBSA)

M507

### Internal Use:

Used for selective isolation of enteric bacteria on the basis of lactose fermentation.

### Composition\*\*

Ingredients	g/ L
Peptone	20.000
Sodium taurocholate	5.000
HM peptone B #	5.000
Sodium chloride	5.000
Lactose	20.000
Litmus	0.500
Agar	15.000
Final pH ( at 25°C)	7.4±0.2

\*\*Formula adjusted, standardized to suit performance parameters

# Equivalent to Beef extract

### Directions

Suspend 70.5 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

Numerous plating media are in use today for the differentiation of lactose-fermenters and lactose non-fermenters. Some of these are selective, whereas others are differential. Some lactose fermenting, gram-negative enteric bacteria can tolerate the inhibitory substances present in the media. These bacteria can be recognized readily by their appearance on selective plates. Litmus Lactose Bile Salt Agar is a modification of Litmus Lactose Agar formulated by Wurtz (1) and is used for the isolation of enteric bacteria. It can be successfully used in place of MacConkey Agar.

LLBSA Medium contains sodium taurocholate, which inhibits the growth of gram-positive microorganisms. Lactose is the fermentable sugar utilized by coliform enteric bacteria leading to production of acid. Peptone and HM peptone B supply the essential nutrients like nitrogen compounds for the growth of enteric bacteria. Sodium chloride maintains the osmotic balance of the medium.

### Type of specimen

Clinical samples- urine, faeces, etc. Water samples

### Specimen Collection and Handling:

For clinical samples, follow appropriate techniques for sample collection, processing as per guidelines.(2,3)

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards.(4)

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro Diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling specimens. Safety guidelines may be referred in individual safety data sheets

### Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Further biochemical and serological test must be performed for confirmation.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Light purple to greyish yellow homogeneous free flowing may contain minute to small particles.

### Gelling

Firm, comparable with 1.5% Agar gel.

### Colour and Clarity of prepared medium

Light purple coloured slightly opalescent gel forms in Petri plates, may have black particles

### Reaction

Reaction of 7.05 w/v aqueous solution at 25°C. pH : 7.4±0.2

### pH

7.20-7.60

### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant	≥50%	red
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	≥50%	deep blue to violet
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	good-luxuriant	≥50%	
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50-100	none-poor	≤10%	
<i>Proteus mirabilis</i> ATCC 25933	50-100	good-luxuriant (no swarming)	≥50%	blue to violet
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	none-poor	≤10%	

Key : \*Corresponding WDCM numbers.

## Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (2,3).

## Reference

1. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
2. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition
3. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
4. Wurtz, 1897, Technique Bacteriologique Paris, Masson.

Revision : 03 /2024



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
[www.cepartner4u.eu](http://www.cepartner4u.eu)



**In vitro diagnostic  
medical device**



**CE Marking**

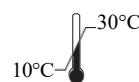

**Storage temperature**



**Do not use if  
package is damaged**

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.