



Moeller Decarboxylase Broth w/ Arginine HCl

M689

Intended Use:

Recommended for differentiation of bacteria on the basis of their ability to decarboxylate L-Arginine hydrochloride.

Composition**

Ingredients	g/ L
Peptone	5.000
HM peptone B #	5.000
Dextrose (Glucose)	0.500
Bromocresol purple	0.010
Cresol red	0.005
Pyridoxal	0.005
L-Arginine hydrochloride	10.000
Final pH (at 25°C)	6.0±0.2

**Formula adjusted, standardized to suit performance parameters

Equivalent to Beef extract

Directions

Suspend 20.52 grams in 1000 ml purified/distilled water. Heat if necessary to dissolve the medium completely. Dispense 5 ml amounts into screw-capped tubes. Sterilize by autoclaving at 15 lbs pressure (121°C) for 10 minutes. Cool the tubed medium in an upright position. Inoculate the tubes and overlay with 2-3 ml of sterile mineral oil.

Principle And Interpretation

Moeller Decarboxylase Broth with Arginine hydrochloride is used for differentiating gram-negative enteric bacilli on the basis of their ability to decarboxylate L-Arginine hydrochloride. Decarboxylase Broth was introduced by Moeller for detecting the production of lysine and ornithine decarboxylase and arginine dihydrolase (1). Prior to Moellers work, bacterial amino acid decarboxylases were studied by Gale (2) and Gale and Epps (3). Decarboxylase media are also recommended by standard methods for identification of bacteria (4,5,6,7).

This medium contains HM peptone B and peptone which provide nitrogenous nutrients for the growth of bacteria. Dextrose is the fermentable carbohydrate and pyridoxal is the co-factor for the decarboxylase enzyme. Bromo cresol purple and cresol red are the pH indicators in this medium. When the medium is inoculated with dextrose fermenting bacteria, the pH is lowered due to acid production which changes the colour of the indicator from purple to yellow.

Acid produced stimulates decarboxylase enzyme. Arginine is first hydrolyzed to ornithine which is then decarboxylated to form putrescine. Formation of the amine putrescine increases the pH of the medium, changing the colour of the indicator from yellow to purple. If the organisms do not produce the appropriate enzyme, the medium remains acidic, yellow in colour. Each isolate to be tested should also be inoculated into the basal medium tube lacking the amino acid. After incubation, a decarboxylase test may show two layers of different colours, yellow and purple. Shake the tube gently before interpreting the results (8).

Inoculated tubes must be protected from air with a layer of sterile mineral oil. Exposure to air may cause alkalization at the surface of the medium which makes the test invalid.

Type of specimen

Isolated Microorganism from clinical and non-clinical samples.

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,9).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling specimens. Safety guidelines may be referred in individual safety data sheets.

Please refer disclaimer Overleaf.

Limitations :

1. Inoculated tubes must be protected from air with a layer of sterile mineral oil. Exposure to air may cause alkalinization at the surface of the medium which makes the test invalid.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to greenish yellow homogeneous free flowing powder

Colour and Clarity of prepared medium

Purple coloured clear solution without any precipitate

Reaction

Reaction of 2.05% w/v aqueous solution at 25°C. pH : 6.0±0.2

pH

5.80-6.20

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for upto 4 days (Inoculated tubes are overlaid with sterile mineral oil).

Organism	Arginine decarboxylation
<i>Citrobacter freundii</i> ATCC 8090	variable reaction
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	negative reaction, yellow colour
<i>Escherichia coli</i> ATCC 25922 (00013*)	variable reaction
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	negative reaction, yellow colour
<i>Proteus mirabilis</i> ATCC 25933	negative reaction, yellow colour
\$ <i>Proteus hauseri</i> ATCC 13315	negative reaction, yellow colour
^ <i>Pseudomonas paraeruginosa</i> ATCC 9027 (00026*)	positive reaction, purple colour
<i>Salmonella</i> Paratyphi A ATCC 9150	delayed positive reaction/positive reaction, purple colour
<i>Salmonella</i> Typhi ATCC 6539	delayed positive reaction /negative reaction, yellow colour
<i>Serratia marcescens</i> ATCC 8100	negative reaction, yellow colour
<i>Shigella dysenteriae</i> ATCC 13313	delayed positive reaction/negative reaction, yellow colour
<i>Shigella flexneri</i> ATCC 12022 (00126*)	delayed positive reaction/negative reaction, yellow colour

Shigella sonnei ATCC
25931

variable reaction

Key : *Corresponding WDCM numbers.

\$ Formerly known as *Proteus vulgaris* ^ Formerly known as *Pseudomonas aeruginosa*

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 15-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,9).

Reference

1. Moeller V., 1955, Acta Pathol. Microbiol. Scand. 36:158.
2. Gale G. F., 1940, Biochem. J., 34:392.
3. Gale and Epps, 1943, Nature, 152:327.
4. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
5. FDA Bacteriological Analytical Manual, 8th Ed., AOAC International, Gaithersburg, Md.
6. Isenberg (Ed.), 1992, Clinical Microbiology Procedures Handbook, Vol. I, ASM, Washington, D. C.
7. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of foods, 5th Ed., American Public Health Association, Washington, D.C.
8. MacFaddin J. F., 2000, Biochemical tests for Identification of Medical Bacteria, 3rd Ed., Lippincott, Williams and Wilkins, Baltimore.
9. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 03/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

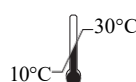

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.