



Nutrient Agar Plate(150mm plate)

MP001CL

Intended use

Recommended a general purpose medium for the cultivation of less fastidious microorganisms.

Composition**

Ingredients	g/ L
Peptone	5.000
Sodium chloride	5.000
HM peptone B#	1.500
Yeast extract	1.500
Agar	15.000
Final pH (at 25°C)	7.4±0.2

**Formula adjusted, standardized to suit performance parameters

- Equivalent to Beef extract

Directions

1. This product is available in a stack of 5 plates.
2. Plan the experiments so as to use all plates at a time.
3. Before taking the plates into the laboratory area where sample is to be processed, follow steps as below.
4. Product has sterile plate containing the medium with three outer protective layers.
5. Place the plates in inverted position while use. Ensure that the cover lids are intact.
6. The outer sterile packaging should be removed carefully by assistant. Personnel with sterile gloves should not touch the outer surface of sterile packaging. Outer packaging must not be placed on sterile surfaces.
7. The product now has two more inner sterile layers. This inner SBS must be handled by person wearing appropriate personal protective equipment (PPE) and sterile gloves. It can be placed on sterile surfaces and opened in sterile or aseptic conditions while working.
8. Remove the second layer and innermost layer and discard.
9. Sterile plate with medium is now ready for use.
10. Aseptically open the plate for inoculation or streaking under aseptic conditions, laminar air flow (LAF) only. Open the plate by holding the lid with three fingers of one hand and base onto the palm as per standard microbiological practice.
11. On completion of inoculation or streaking close the lid carefully and place the petriplate in inverted position.
12. Incubate as desired.
13. Follow good lab practices for procedures and disposal.

Principle And Interpretation

Nutrient media are basic culture media used for maintaining microorganisms, cultivating fastidious organisms by enriching with serum or blood and are also used for purity checking prior to biochemical or serological testing (1,2). Nutrient Agar is ideal for demonstration and teaching purposes where a more prolonged survival of cultures at ambient temperature is often required without risk of overgrowth that can occur with more nutritious substrate. This relatively simple formula has been retained and is still widely used in the microbiological examination of variety of materials and is also recommended by standard methods. It is one of the several non-selective media useful in routine cultivation of microorganisms (3,4). It can be used for the cultivation and enumeration of bacteria which are not particularly fastidious. Addition of different biological fluids such as horse or sheep blood, serum, egg yolk etc. makes it suitable for the cultivation of related fastidious organisms. Peptone, HM peptone B and yeast extract provide the necessary nitrogen compounds, carbon, vitamins and also some trace ingredients necessary for the growth of bacteria. Sodium chloride maintains the osmotic equilibrium of the medium.

Type of specimen

Clinical samples - faeces, urine ; Food and dairy samples; Water samples

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6). For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (3,4,7). For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (8). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's requirement.
3. Further serological and biochemical testing is required for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile Nutrient Agar in 150 mm disposable plates with smooth surface and absence of black particles/cracks/bubbles **Colour of medium**

Light amber coloured medium

Quantity of medium

85 ml of medium in 150 mm disposable plates.

pH

7.20-7.60

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant	≥70%
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	good-luxuriant	≥70%
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	≥70%
<i>Staphylococcus aureus</i> subsp.aureus ATCC 25923 (00034*)	50-100	good-luxuriant	≥70%
<i>Streptococcus pyogenes</i> ATCC 19615	50-100	good-luxuriant	≥70%
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	good-luxuriant	≥70%

<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	good-luxuriant	$\geq 70\%$
<i>Yersinia enterocolitica</i> ATCC 9610 (00038*)	50-100	good-luxuriant	$\geq 70\%$
<i>Yersinia enterocolitica</i> ATCC 23715 (00160*)	50-100	good-luxuriant	$\geq 70\%$

Key : (*) Corresponding WDCM numbers.

Storage and Shelf Life

On receipt store between 20-30°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

Reference

- Lapage S., Shelton J. and Mitchell T., 1970, Methods in Microbiology', Norris J. and Ribbons D., (Eds.), Vol. 3A, Academic Press, London.
- MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Lippincott, Williams and Wilkins, Baltimore.
- American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
- Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
- Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
- Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
- Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.

Revision : 00/2024

Packaging:

MP001CL - 20PT -Mueller Hinton Agar Plate (150mmX 5Plates)


MP001CL - 50PT -Mueller Hinton Agar Plate (150mmX 5Plates)




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



Triple sterile
barrier system with
protective packaging
outside

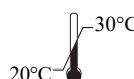

Storage temperature

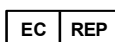

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



In vitro diagnostic
medical device



Do not use if
package is damaged



Do not re-use



CE Marking



Contains biological
material of animal origin

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.