

Sabouraud Chloramphenicol Agar Plate

MP1067

Intended use

Recommended for the selective cultivation of yeasts and moulds from clinical and non-clinical samples.

Composition**

Ingredients	g / L
Tryptone	5.000
Peptone	5.000
Dextrose (Glucose)	40.000
Chloramphenicol	0.050
Agar	15.000
Final pH (at 25°C)	5.6±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

Caution: Some pathogenic fungi may produce infective spores which are easily dispersed in air, so examination should be carried out in safety cabinet.

Principle And Interpretation

Sabouraud Chloramphenicol Agar is cited as Medium C and recommended for cultivation of yeasts and moulds. This medium was described originally by Sabouraud (1) for the cultivation of fungi, particularly useful for the fungi associated with skin infections. The medium is often used with antibiotics such as Chloramphenicol (2) for the isolation of pathogenic fungi from materials containing large numbers of fungi or bacteria.

Tryptone and peptone provide nitrogenous and carbonaceous compounds, long chain amino acids, and other essential growth nutrients. Dextrose provides an energy source. Chloramphenicol inhibits a wide range of Gram-positive and Gram-negative bacteria which makes the medium selective for fungi (3). The low pH favors fungal growth and inhibits contaminating bacteria from clinical specimens (4).

Type of specimen

Clinical samples - skin scrapings, nail scrapings, etc. ; Food and dairy samples.

Specimen Collection and Handling:

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (5,6,7). For clinical samples follow appropriate techniques for handling specimens as per established guidelines (8,9). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use. For professional use only. Read the label before opening the pack. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

- 1.Certain pathogenic fungi may show poor growth on this medium.
- 2.Presence of chloramphenicol may inhibit certain pathogenic fungi.
- 3.Overheating of the medium may result in low productivity and softening of gel.
- 4.Individual organisms differ in their growth requirement and may show variable growth patterns on the medium
- 5.Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
- 6.It is recommended to store the plates at 24-30°C to avoid condensation.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile Sabouraud Chloramphenicol Agar in 90 mm disposable plates with smooth surface and absence of black particles/cracks/bubbles

Colour of medium

Light amber coloured medium

Quantity of medium

25 ml of medium in 90 mm disposable plates.

pH

5.40-5.80

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after an incubation at 20-25°C for 48-72 hours (Incubate for 7 days for *Trichophyton* species).

Organism	Inoculum (CFU)	Growth	Recovery
# <i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	50-100	good-luxuriant	
<i>Candida albicans</i> ATCC 10231 (00054*)	50-100	good-luxuriant	≥50%
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 ³	inhibited	0%
\$ <i>Lactobacillus paracasei</i> ATCC 334	≥10 ³	inhibited	0%
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	50-100	good-luxuriant	≥50%
<i>Trichophyton rubrum</i> ATCC 28191	50-100	good-luxuriant	
<i>Escherichia coli</i> ATCC 8739 (00012*)	≥10 ³	inhibited	0%

Key : (*)Corresponding WDCM numbers. (\$) Formerly known as *Lactobacillus paracasei*
 (#)Formerly known as *Aspergillus niger*

Storage and Shelf Life

On receipt store between 20-30°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (8,9).

Reference

- Sabouraud K., 1892, Ann. Dermatol. Syphilol, 3:1061.
- Ajello L., 1957, J. Chron. Dis., 5:545.
- Lorian (Ed.),1980, Antibiotics In Laboratory Medicine, Williams and Wilkins, Baltimore.
- Murray, P. R 2005, In Manual of Clinical Microbiology, 7th ed., ASM, Washington, D.C.
- Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.

6. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
7. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
8. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
9. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 04/2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



IVD *In vitro* diagnostic
medical device

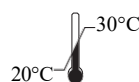

Storage temperature

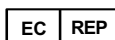

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.