

SS Agar (Salmonella Shigella Agar) Plate

MP108

Intended Use:

Recommended for differential and selective isolation of *Salmonella* and *Shigella* species from pathological specimens, suspected foodstuffs etc.

Composition**

Ingredients	g / L
Peptone	5.000
HM peptone B #	5.000
Lactose	10.000
Bile salts mixture	8.500
Sodium citrate	10.000
Sodium thiosulphate	8.500
Ferric citrate	1.000
Brilliant green	0.00033
Neutral red	0.025
Agar	15.000
Final pH (at 25°C)	7.0±0.2

**Formula adjusted, standardized to suit performance parameters

- Equivalent to Beef extract

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

Principle And Interpretation

SS Agar medium is recommended as differential and selective medium for the isolation of *Salmonella* and *Shigella* species from pathological specimens (1) and suspected foodstuffs (2,3,4,5) and for microbial limit test (6). SS Agar is a moderately selective medium in which gram-positive bacteria are inhibited by bile salts, brilliant green and sodium citrate. Peptone, HM peptone B provides nitrogen and carbon source, long chain amino acids, vitamins and essential growth nutrients. Lactose is the fermentable carbohydrate. Brilliant green, bile salts and thiosulphate selectively inhibit gram-positive and coliform organisms. Sodium thiosulphate is reduced by certain species of enteric organisms to sulphite and H₂S gas and this reductive enzyme process is attributed by thiosulphate reductase. Production of H₂S gas is detected as an insoluble black precipitate of ferrous sulphide, formed upon reaction of H₂S with ferric ions or ferric citrate, indicated in the center of the colonies.

The high selectivity of Salmonella Shigella Agar allows the use of large inocula directly from faeces, rectal swabs or other materials suspected of containing pathogenic enteric bacilli. On fermentation of lactose by few lactose-fermenting normal intestinal flora, acid is produced which is indicated by change of colour from yellow to red by the pH indicator-neutral red. Thus these organisms grow as red pigmented colonies. Lactose non-fermenting organisms grow as translucent colourless colonies with or without black centers. Growth of *Salmonella* species appears as colourless colonies with black centers resulting from H₂S production. *Shigella* species also grow as colourless colonies which do not produce H₂S.

Type of specimen

Clinical: faeces, rectal swabs, etc.; Suspected food stuffs.

Specimen Collection and Handling

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (2,3,4). For clinical samples follow appropriate techniques for handling specimens as per established guidelines (7,8). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

- 1.The medium is highly selective and may be toxic to certain *Salmonella* or *Shigella* species. Hence it is recommended to use to inoculate plates of less inhibitory media parallel to SS Agar, such as Hektoen Enteric Agar (M467) or Deoxycholate Citrate Agar (M065) for easier isolation of *Shigella* species (3).
- 2.Individual organisms differ in their growth requirement and may show variable growth patterns on the medium
- 3.Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile SS Agar (Salmonella Shigella Agar) in 90 mm disposable plates with smooth surface and absence of black particles/cracks/bubbles

Colour of medium

Reddish orange coloured medium

Quantity of medium

25 ml of medium in 90 mm disposable plates.

pH

6.80-7.20

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	fair	20-30%	cream pink
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	fair	20-30%	pink with bile precipitate
<i>Salmonella Choleraesuis</i> ATCC 12011	50-100	good-luxuriant	>=50%	colourless with black centre
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	>=50%	colourless with black centre
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	none-poor	<=10%	colourless
<i>Proteus mirabilis</i> ATCC 25933	50-100	fair-good	30-40%	colourless, may have black centre
<i>Shigella flexneri</i> ATCC 12022 (00126*)	50-100	good	40-50%	colourless
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	good-luxuriant	>=50%	colourless with black centre
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	good-luxuriant	>=50%	colourless with black centre

Key : (*) Corresponding WDCM numbers. (#) Formerly known as *Enterobacter aerogenes*

Storage and Shelf Life

On receipt store between 20-30°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (7,8).

Reference

- 1.Lennette and others (Eds.), 1985, Manual of Clinical Microbiology, 4th ed., ASM, Washington, D.C.
- 2.Lipps WC, Braun-Howland EB, Baxter TE,eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
- 3.Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- 4.Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed.,APHA Inc., Washington, D.C.
- 5.Williams S., (Ed.), 2005, Official Methods of Analysis of the Association of Official Analytical Chemists, 19th Ed.,AOAC,Washington, D.C.
- 6.The United States Pharmacopoeia, 2006, USP29/NF24, The United States Pharmacopoeial Convention. Rockville, MD.
- 7.Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
- 8.Jorgensen,J.H., Pfaller , M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision :03/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



IVD *In vitro diagnostic
medical device*

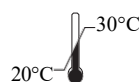

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.