

HiCrome™ Salmonella Agar Plate

MP1296

Intended Use:

Recommended for the isolation and differentiation of *Salmonella* species from coliforms by chromogenic method.

Composition**

Ingredients	g / L
Peptone	6.000
Yeast extract	2.500
Bile salts mixture	1.000
Chromogenic mixture	5.400
Agar	13.000
Final pH (at 25°C)	7.7±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

Principle And Interpretation

Salmonella species have been isolated from humans and almost all animals throughout the world. They cause many types of infections from mild, self-limiting gastroenteritis to life threatening typhoid fever. *Salmonella* Typhi and *Salmonella* Paratyphi A & B cause gastroenteritis, bacteremia and enteric fever, *Salmonella* Choleraesuis causes gastroenteritis and enteric fever, especially in children. *Salmonella* Typhimurium is the most frequently isolated serotype of *Salmonella*. *Salmonella* is a cause of food poisoning (1).

HiCrome™ Salmonella Agar is a modification of the original formulation of Rambach (2) and is used for the differentiation of *Salmonella* species from other enteric bacteria. Rambach formulation differentiates *Salmonella* based on propylene glycol utilization and presence of a chromogenic indicator. However, HiCrome™ Salmonella Agar medium uses only a chromogenic mixture for identification and differentiation of *Salmonella* species.

Peptone and yeast extract provides nitrogenous, carbonaceous compounds and other essential growth nutrients. *Escherichia coli* and *Salmonella* are easily distinguishable due to their colony characteristics. *Salmonella* forms light purple coloured colonies with a purple halo.

E.coli and other β-glucuronidase positive organism exhibits a characteristic blue colour, due to presence of the enzyme β-glucuronidase. Other organisms form colourless colonies. The characteristic light purple and blue colour is due to the chromogenic mixture (3). Bile salts mixture inhibits gram-positive organisms.

Conventional method employees the H2S production property for *Salmonella* detection which is also exhibited by other non *Salmonella* species such as *Citrobacter*, *Proteus*, etc. Hence further biochemical confirmation is required for further identification. *Salmonella* species isolated from food or clinical samples exhibit light purple colour with halo due to the specific enzyme substrate reaction.

Type of specimen

Clinical samples: faeces, urine, etc.; Food Samples; Water samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (6). For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (7). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use. For professional use only. Read the label before opening the pack. Wear protective gloves/ protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets. Precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. The medium is selective for *Salmonella* may not support the growth of other microorganisms.
2. Most of the *Salmonella* strains show light purple colonies except few which may show colorless colonies.
3. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium
4. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
5. Final confirmation of suspected colonies must be carried out by serological and biochemical tests.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile HiCrome Salmonella Agar in 90mm disposable plate with smooth surface and absence of black particles/cracks/bubbles.

Colour of medium

Cream coloured opalescent medium

Quantity of medium

25 ml of medium in 90 mm disposable plate

pH

7.50-7.90

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of Colony
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	$\geq 10^3$	inhibited	0%	
** <i>Bacillus spizizenii</i> ATCC 6633 (00003*)	$\geq 10^3$	inhibited	0%	
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	$\geq 50\%$	blue
## <i>Proteus hauseri</i> ATCC 13315	50-100	good	40-50%	colourless
<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50-100	luxuriant	$\geq 50\%$	light purple w/ halo
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	luxuriant	$\geq 50\%$	light purple w/ halo
<i>Salmonella</i> Typhi ATCC 6539	50-100	good-luxuriant	$\geq 50\%$	light purple w/ halo

Key- (*) Corresponding WDCM numbers

(**) Formerly known as *Bacillus subtilis* subsp. *spizizenii*

(##) Formerly known as *Proteus vulgaris*

Storage and Shelf Life

On receipt store between 2-8°C. Use before expiry date on the label. Product performance is best if used within stated expiry period. .

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference

- 1.Murray P. R., Baron J. H., Pfaller M. A., Jorgensen J. H. and Tenover F. C., (Ed.), 2003, Manual of Clinical Microbiology, 8th Ed., American Society for Microbiology, Washington, D.C.
- 2.Rambach A., 1990, Appl. Environ. Microbiol., 56:301.
- 3.Greenwald R., Henderson R. W. and Yappan S., 1991, J. Clin. Microbiol., 29:2354.
- 4.Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
- 5.Jorgensen, J.H., Pfaller , M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
- 6.Salfinger Y., and Tortorello M.L. , 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
- 7.Lipps WC, Braun-Howland EB, Baxter TE,eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.


Revision : 02/ 2024




HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



**In vitro diagnostic
medical device**

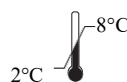

Storage temperature

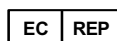

CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.